

i6000™ Infinity Automated Staining System

Multifunctional Staining System for Research

OPERATOR'S MANUAL

i6000™ Infinity Automated Staining System (Software Version 4.3r)
Doc. No. 953-0053.0, Rev. D, ECO: 303-SEZ

CE

For Research Use Only



48810 Kato Road,
Suite 100E & 200E,
Fremont, CA 94538
Tel : +1 (800) 421-4149
Fax : +1 (510) 824-1490
E-mail: support@biogenex.com

 Plot No: 4, Survey No: 656/A,
Aerospace and precision engineering industries SEZ.
Adibatla (V), Ibrahimpatnam (M), Ranga Reddy (D)
Telangana, India.
Pin: 501510.

Technical Support:

For USA and Canada:

Tel: 1.800.421.4149

Fax: 1.925.866.2594

For outside USA and Canada

Email: support@biogenex.com

www.biogenex.com

i6000™ Infinity Automated Staining System

Copyright ©2010 by BioGenex Laboratories Inc. All Rights Reserved



Emergo Europe

Prinsessegracht 20

2514 AP The Hague

The Netherlands

Tel: (+31) 70-345-8570

Fax: (+31) 70-346-7299

Table of Contents

Section No	Title	Page number
1.0	Introduction to i6000 InfinityTM Automated Staining System	6
1.1	The Instrument	6
1.1.1	Intended Use	6
1.1.2	Automatic Universal Staining	6
1.2	General Principles of Operation	6
1.3	Precautions and Warning	6
2.0	Instrument Overview	9
2.1	Accessory Kit	10
2.2	Instrument Setup	11
2.2.1	Loading the Slides	12
2.2.2	Loading Reagents	13
3.0	Getting Started	14
3.1	Logging On	14
4.0	Tool Bars and Menu	18
4.1	Tool Bar	18
4.2	Global Settings	20
4.2.1	Buffer Options	20
4.2.2	Liquid Sensor Option	20
4.2.3	Use Polymer Detection Kit for IHC	20
4.2.4	Stain Pattern	20
5.0	Executing a Run	22
5.1	Inserting Protocol	22
5.2	Copying Protocols	23
5.3	Edit Slide Parameters	24
5.4	Deleting Slides	25
5.5	Starting the Run	25
5.6	Reagent Screen	27
5.7	Delayed run	28
5.8	Slide Information	30
5.9	Slide Priority	30
5.10	Run Information	31
5.11	Runtime Properties	31
5.12	Inserting Continuous Slides	32
5.13	Inserting STAT Slides	33
5.14	Run Completion	35
6.0	Protocols	36
6.1	Creating a New Protocol	36
6.2	Copying a Protocol	39
6.3	Edit Protocol	40

6.4	Revert Protocol	41
6.5	Deleting a Protocol	42
7.0	Reagents	44
7.1	Creating a New Reagent	44
7.2	Editing a Reagent Description	44
8.0	Global Settings and Other Tools	45
8.1	User Manager	45
8.2	Edit User Profiles	45
8.3	Delete User Profile	46
8.4	Change Password	46
8.5	Physician Names	46
8.6	Specimen Types	47
8.7	Color Scheme	47
9.0	Pause or Stop Function	48
9.1	Emergency Stop	48
9.2	Software Pause	48
9.3	Software Stop	49
10.0	Reports	50
10.1	Starting the Reports Module	50
10.2	Export a Report	52
11.0	Calibration and Diagnostics	53
11.1	One Point Calibration	54
11.2	Multi Point Calibration	55
11.3	Left Pipette Rack	55
11.4	Right Pipette Rack	56
11.5	Left Reagent Rack	57
11.6	Right Reagent Rack	58
11.7	Slide Rack	58
11.8	Prime Position	59
11.9	Pipette Tip Disposal	59
11.10	Vial Bottom Calibration	59
11.11	Diagnostics	61
11.12	Position	62
11.13	Operations	62
11.14	Sensors and Valves	63
12.0	Error Conditions	63
12.1	System Door Error	64
12.2	Pipette Errors	64
12.2.1	Pipette Tip Errors at Box	64
12.2.2	Pipette Tip Errors at Reagent	65
12.2.3	Pipette Tip Errors at Slide	66
12.2.4	Pipette Tip Errors post Dispense	67
12.3	Reagent Errors	68
12.3.1	Insufficient Reagent Vials Error	69

12.4	System Pressurization Errors	69
12.5	Substrate and Label Error	69
Appendix-1	System Specifications	71
Appendix-2	Service Plan	72
Appendix-3	Software License	73
Appendix-4	Software License Agreement	75
Appendix-5	Daily/ Monthly Cleaning	78
Appendix-6	Instrument Packaging Procedure for i6000 Diagnostics/ Infinity Automated Staining System	80
Appendix-7	Pre-Installation Checklist for i6000™ Systems	84
Appendix-8	Instrument Installation Instructions for i6000 Diagnostics/ Infinity Automated Staining Systems	89

1.0 Introduction

i6000™ Infinity - Automated Staining System

1.1 The Instrument

1.1.1 Intended Use

The i6000™ Infinity is intended for research use in Immunohistochemistry (IHC), In-Situ Hybridization (ISH), Special Stains and other histological and histochemical staining procedures. The results must be interpreted by a qualified pathologist.

1.1.2 Automated Universal Staining

The i6000™ Infinity is a complete and reliable system used for histochemical, immunohistochemical, and immunocytochemical staining in addition to certain In situ hybridization protocols. This fully automated, walk-away system provides excellent and consistent staining results between assays and can stain up to up to 60 slides in one single run. The quick and efficient robotic execution of protocols saves time, money, and labor thereby allowing laboratory workers to perform other tasks. The i6000™ Infinity combines BioGenex's superior primary antibodies, reagents and detection systems with a cutting-edge system automation and the user-friendly i6000™ Windows XP/7 Software. BioGenex provides a complete set of reagents, buffers, hardware, and software, optimized for inter-assay consistency and accuracy of staining. The i6000™ Infinity operates in an open "user defined" mode. The user defined mode allows users to develop "home-brewed" protocols for their unique applications, thus meeting the changing needs of reference and clinical laboratories, and the basic user defined community.

1.2 General Principles of Operation

The i6000™ Infinity consists of an XYZ robotic arm with a unique staining head that picks-up and discards standard disposable pipette tips, and precisely applies reagents to tissue sections thereby eliminating the possibility of cross-contamination. The integrated wash and blow head rinses the slides after each reagent-incubation and blows off the excess buffer from the slides at the end of the rinsing step. It is designed to produce maximal spreading of reagents on tissue sections when used in conjunction with BioGenex immunostaining reagents. The i6000™ Infinity is designed to complete a staining run of 60 slides in less than three and half-hours.

1.3 Precautions and Warning

The i6000™ Infinity must be unpacked, inspected, and installed by an authorized BioGenex representative. The serial number should be recorded and the accessories should be checked against the accessory list. The i6000™ Infinity presents no hazards to operators if

operated according to the instructions in this manual. Read the following safety precautions before applying power to the instrument.

The i6000™ Infinity may be used only for research uses specified in this manual and for appropriate user defined purposes. Any other use is prohibited. Clinical interpretation of any staining results from the i6000™ Infinity operation is solely the responsibility of the user.

The clinical interpretation of any positive signal or its absence should be evaluated within the context of clinical presentation, morphology and other histopathological criteria. The clinical interpretation of any positive signal or its absence should be complemented by morphological studies using proper positive and negative internal and external controls, as well as other diagnostic tests.

Do not operate the i6000™ Infinity in an environment where inflammable vapors may be present. Operating AC-line-powered equipment in such an environment may cause an explosion.



Fuse type: GDB5A, 250V

Warning: For continued protection against fire hazard, replace only with same type and rating of fuse



When replacing any fuse on the i6000™ Automated Staining System, make sure that a fuse of the same type and rating is used. Using the wrong fuse may lead to fire or instrument failure.

The AC-line cord provided with the i6000™ Infinity is intended to ground the chassis to help prevent shock and injury to personnel. The i6000™ Infinity meets all national and local electrical codes that require a three-conductor AC outlet that connects the third wire of the line cord to the earth ground. Do not in any way defeat the grounding connection of this instrument.

Use of the i6000™ Infinity may require the use of hazardous chemicals. Please refer to the reagent manufacturer's instructions and be sure to follow all applicable local regulations for the use, handling, storage, and disposal of any hazardous chemicals and waste. Liquids used on the i6000™ Infinity may present a slip hazard if spilled on the floor. Excessive jarring of the instrument while the lid is open may cause the lid to close unexpectedly.

Test for ground continuity between the system chassis and adjacent equipment or metal plumbing to ensure that the operator cannot become a conductor between power ground and building ground. Never tamper with the 3-pin power plug.

Do not service the system while power is ON. Refer all servicing to a qualified technician.



Accessing the i6000™ Infinity while the robotics are in motion, may cause bodily injury. Ensure that the run is paused prior to accessing the working space.

The i6000™ Infinity is equipped with a door switch which when open will stop the movement of the robotic arm.



This symbol identifies the location of the earth grounding point on the instrument.

Warning Instructions

Power outlet should be near the equipment and it is easily accessible for easy disconnection of power from the equipment in case of fault.

Important safety instructions:

Misuse of electrical components of i6000™ can cause electrocution, burns, fire and other Hazards.


Read this before using the i6000™

- Connection to MAINS supply: Where protective earthing is required, plug the i6000™ into a supply Qualified electrician or authorized service engineer.
- Use the i6000™ only for the purpose described in the instructions for use.
- Do not use accessories which are not supplied or recommended by BioGenex.
- Do not use the i6000™ if it is not working properly, or if it has suffered any damage.
- Do not let the i6000™ or its flexible cord come into contact with surfaces which are too hot to touch.
- Do not block air openings nor place i6000™ on a soft surface which might block them, and keep air openings free from lint, hair, fluff, etc.
- Do not place anything on top of the i6000™.
- Unless specifically instructed to do so by the instructions for use do not drop or put anything into any opening in the i6000™, or into any hose or coupling.
- Do not use the i6000™ out of doors.

For any support, troubleshooting and repair contact BioGenex Technical Support.

2.0 Instrument Overview i6000™Infinity

The BioGenex i6000™ Infinity is a combination of hardware and software integrated with special reagents to automate IHC and special stain methods. This section describes the hardware and software specifications and the accessories provided with the system.

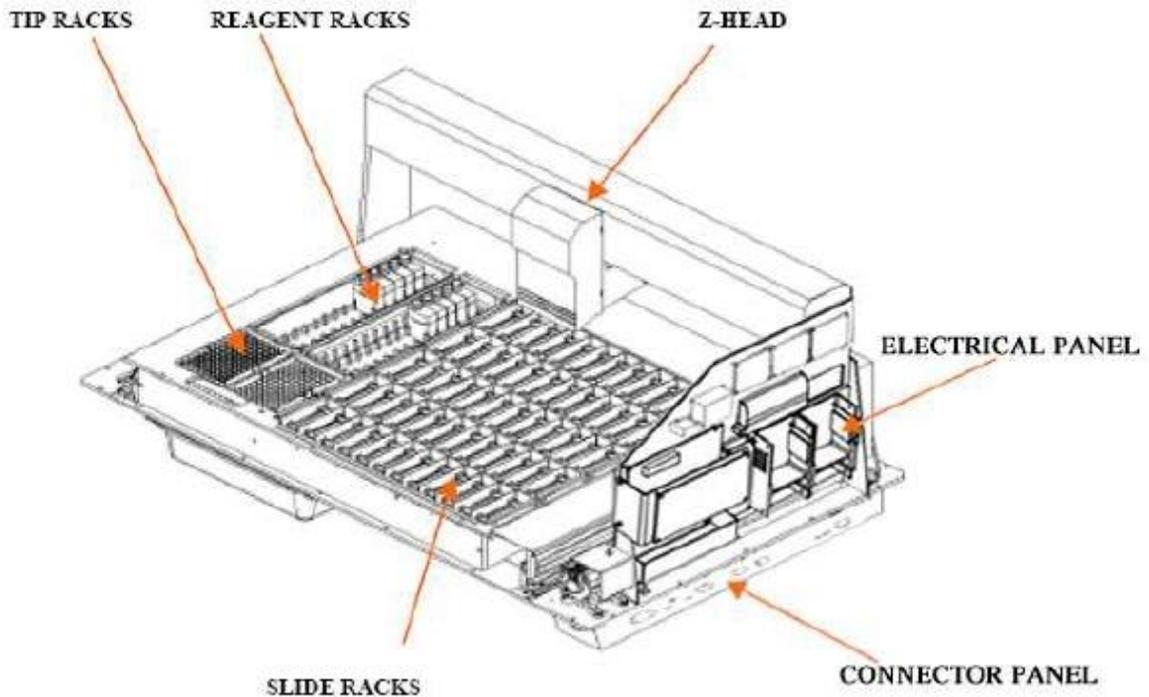

SLIDE RACKS: Consists of 5 racks where each rack can hold up to 12 slides, total 60 slides per run.

TIP RACKS: Consists of 2 racks where each can hold up to 96 tips, total 192 tips in each run. To run the system, use only BioGenex tips.


REAGENT RACKS: Consists of 2 reagent racks, each can hold up to 30 vials, total 60 vials in each run.

Z HEAD: An XYZ robot arm with a unique staining head that picks up and discards standard pipette tips, and precisely applies reagents to tissue sections eliminating any possibility of cross-contamination. The integrated wash and blow head rinses the slides after each reagent incubation and blows the excess buffer from the slide at the end of the rinsing step.

Y-AXIS MECHANISM: This will move the Z Head from front to back positions.

X-AXIS MECHANISM: This will move the Z Head from the left to right positions.

Z-AXIS MECHANISM: This will move the Z Head from the up to down positions.

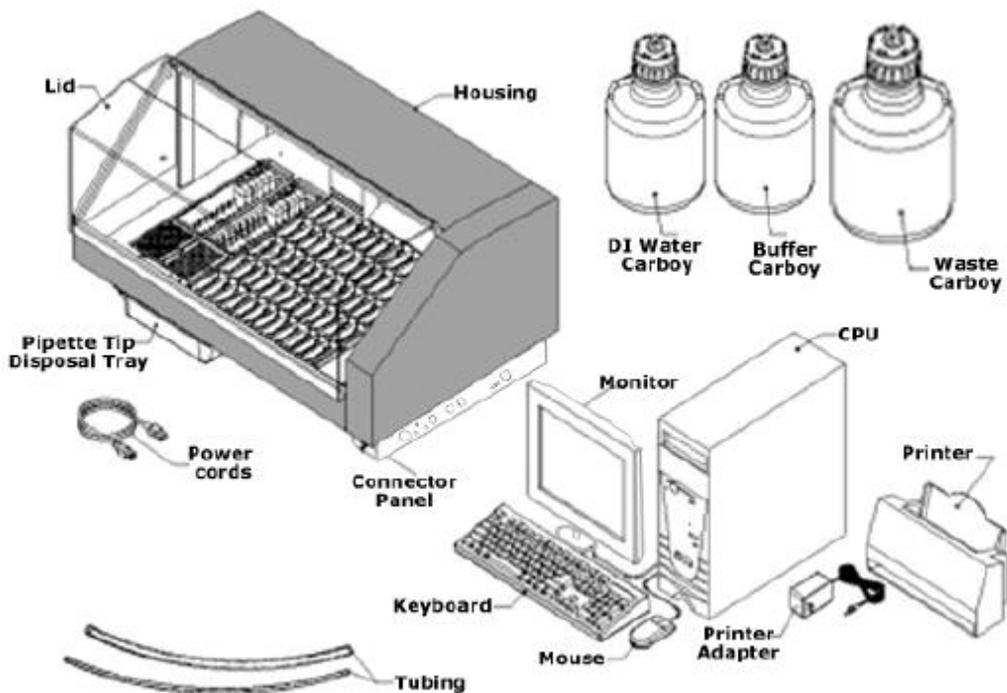

i6000™ Infinity Accessories

2.1 Accessory Kit

BioGenex provides the following accessories along with the i6000™ Automated Staining System.

1.	i6000™ Instrument	01
2.	Carboy and Tubing	
	Air Tubing	01
	Waste Carboy Cap/20L Waste Carboy/Tubing/Label	01
	DI Water Carboy Cap/ Tubing/Label	01
	Buffer Carboy Cap/Tubing/Label	01
	10L Carboy	02
3.	Power Cord	01
4.	Computer Accessories	
	Computer	01
	Monitor	01
	Keyboard	01
	Mouse	01
	Printer	01*
	Printer Cable	01*

(*order separately)


2.2 Instrument Setup

The i6000™ Infinity is installed by authorized BioGenex Service personnel. Once installed and tested on site, the system is intended to remain in that location. If it is necessary to relocate the i6000™ Infinity to another site, contact BioGenex to prepare the system for transport and reinstallation. Failure to properly prepare the system for transport can result in damage to the automated mechanisms.

To use the system, check that the Buffer Carboy is full, and ensure that the tubing from the Buffer Carboy is securely connected to the instrument. An IHC staining run of 60 slides with a factory protocol will use approximately 7.3 liters of buffer.

Note: For best results, be sure to use Super Sensitive Wash Buffer (Cat No. HK583-5K). This buffer has been optimized to maximize the spreading of reagents on the tissue sections.


Make sure the cap O-rings (circular black rubber washers used to prevent air leakage from the carboys) are in place and that the caps are secured tightly to the carboys.

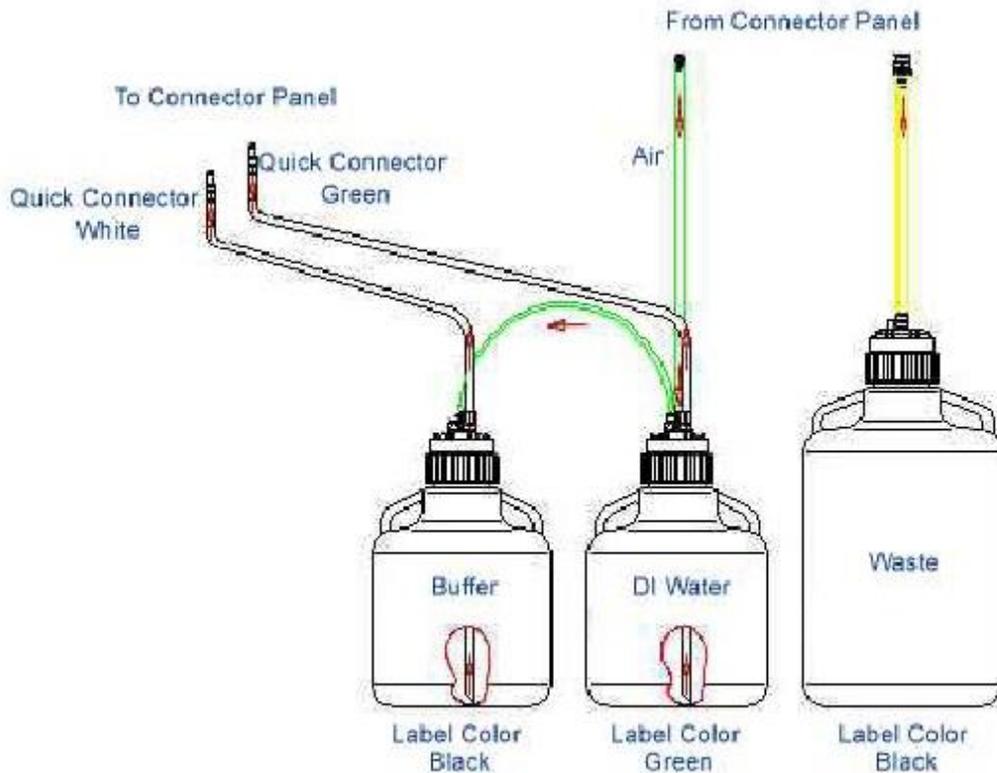

Ensure that the waste carboy is connected to the waste outlet on the side towards the front of the instrument and that the drain line is not kinked.

Make sure the pipette tip disposal tray is in place in the waste tray drawer under the left side of the instrument.


To install a tray of pipette tips in the unit, first place the front edge of the tip tray into position, then press the back side tray down into position. Make sure that the tray is pressed all the way down into the recess. The tip adapter on the Z-head will pick up the first tip at the rear left-hand corner of the tray.

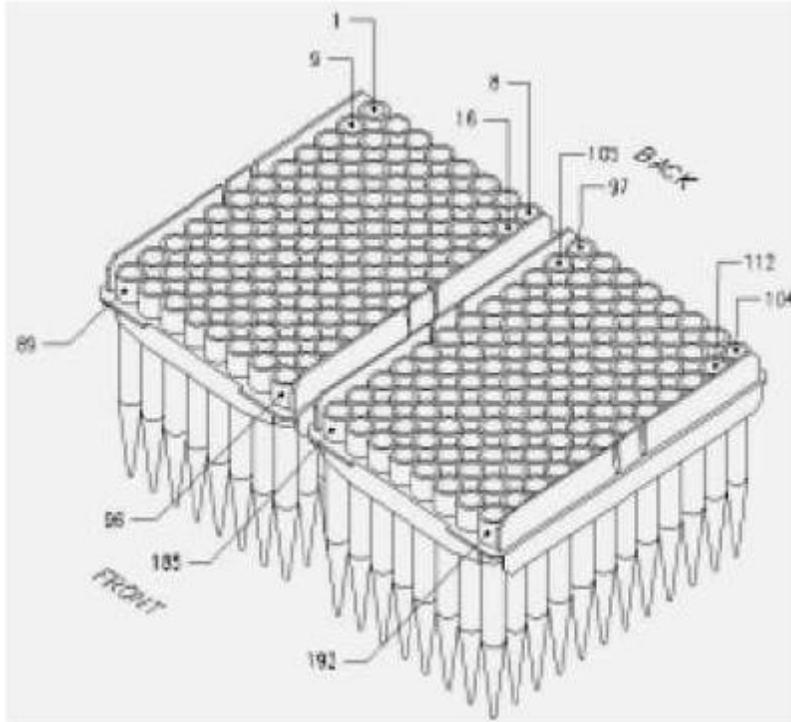

2.2.1 Loading the Slides

Slides can be loaded into the slide racks while they are in the instrument, or the racks can be removed from the System to make loading easier. If you remove the racks, it is recommended that you keep them properly oriented with respect to the instrument, i.e., with the open side towards the front.

For an open run, place the slides into their appropriate positions on the Slide Racks according to your Slide Set-Up Map. Slide Racks have one completely open side for easy loading and one notched side. The frosted end of the slide should be placed at the closed side of the Slide Rack.

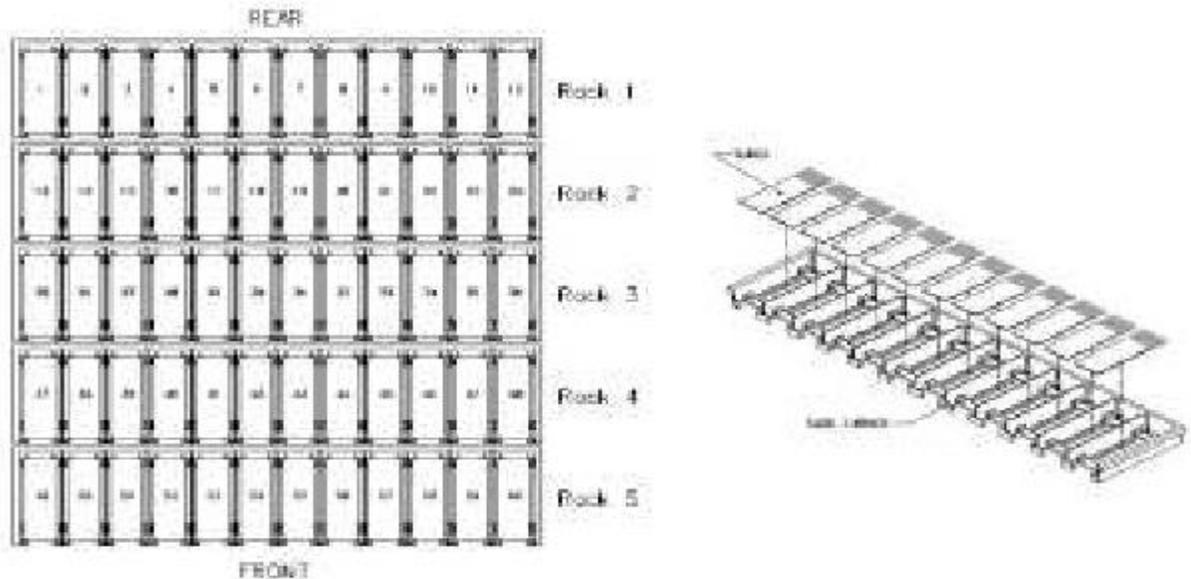

If a rack was removed from the instrument for slide loading, return it to its proper position in the system. The slides must lie flat in the carriers, i.e. the bottom of the slide against the positioning knobs, not on them. When placed in the system, the open side of the Slide Carrier should face the front side of the instrument. Cover the slides with buffer to prevent them from drying out.

2. 2.2 Loading Reagents


The i6000™ Infinity can hold up to 60 different reagents. There are two removable reagent racks. Each can hold 30 vials. The left-hand rack holds vials 1-30 while the right-hand rack holds vials 31-60. Remove the reagent racks from the System to make vial loading easier.

Caution: Use only BioGenex supplied 20ml vials. Improper fluid detection could occur with incorrect vials.

In the diagnostics version, the system will detect the location of the vials. However, orient the vials with the opening facing the left side.

It is very important to place enough reagents into each assigned position and load them in the exact order shown in the reagent map.

When the Reagent Rack is filled, place it in the instrument. Make sure that the rack is oriented properly, i.e. vial position 1 and 31 are installed towards the rear of the instrument. Make sure that the rack is firmly secured in place.

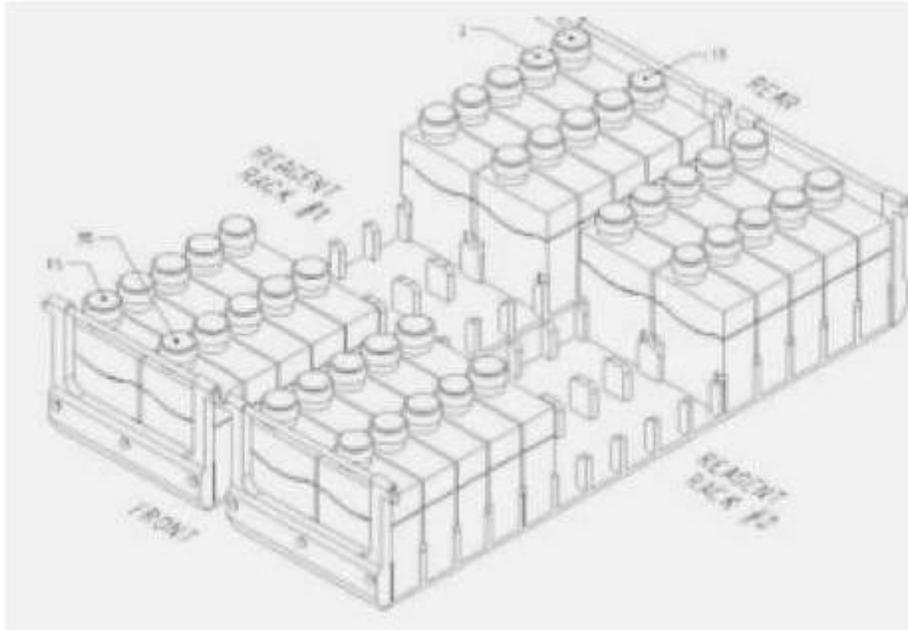

3.0 **Getting Started** **i6000™ Infinity**

3.1 **Logging On**

The i6000™ software runs on the Windows XP®/7 operating system. It was designed to provide the user, a very simple-to-use interface. This section will provide the details of the basic controls used to operate the i6000™ System.

- Switch-on the instrument and the computer.
- 'AutoStainer' is the default user name and enter 'optimax' as password to log on. Password is case sensitive. Use only lower case.
- Make sure that no error messages are displayed as Windows XP/7 starts.
- The Windows® XP/7 desktop will be displayed.
- Double click the i6000™ icon on the desktop.



When the Log on to i6000™ window is displayed, type in the user name and the password. Choose a language of operation from the drop down menu.



The default password for the Administrator is '**BioGenex**'. The password is not case sensitive. Once logged in as the Administrator, additional user profiles can be created. Users other than the administrator will have only limited access as defined by the Administrator.

The <**Welcome to the i6000**> System dialogue box will appear. There are four options: Run i6000™, Print Reports, Calibrate i6000™ and Exit.

Click on:

- **Run i6000™** option to run the instrument and compose runs.
- **Print Reports** to print run reports.
- **Calibrate i6000™** to calibrate the instrument and
- **Exit** to exit the i6000™ application.



Click on the **Run i6000™** button to open the **<Selection of Assay>** dialog box.



Click on the **Run i6000™** button to open the **<Selection of Assay>** dialog box



Click on the **Print Reports** button to open the **<Welcome to Reports i6000>** dialog box




Click on the <Calibrate i6000> button to open the < i6000 Calibration & Diagnostic Tool> dialog box




Click on the **Exit** button to close the application

Select the assay of interest from the three options available – IHC Protocols, ISH Protocols and Special Stains Protocols.

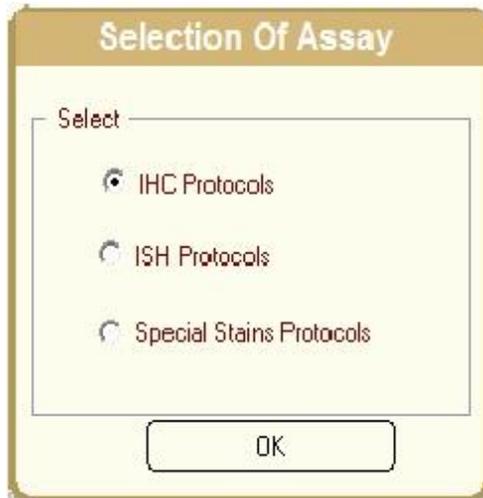

Click 'OK' button to open the <Open Format – i6000 Automated Staining System [Design Window]>

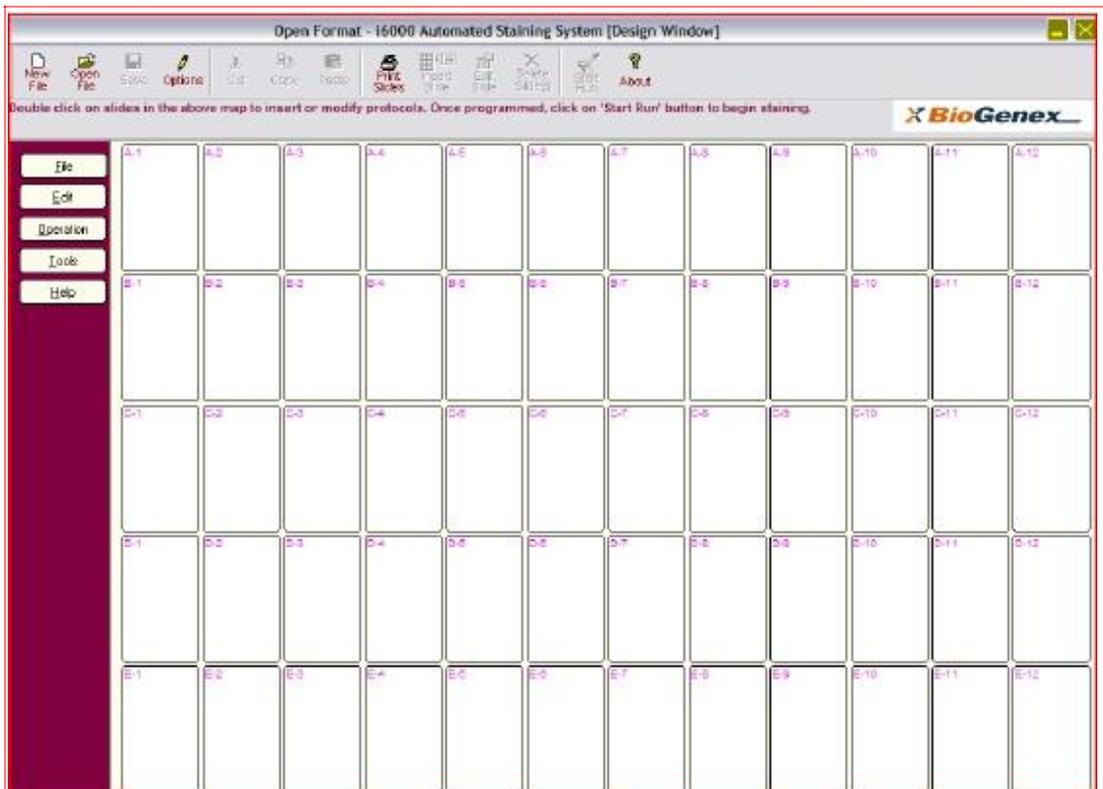

4.0 Tool Bars and Menu

4.1 Tool Bar

Initially when the Design Window is opened, the 'New File', 'Open File', 'Options', 'Print Slides' and 'About' buttons are active.



Once any slide is inserted, the 'Save', 'Cut', 'Copy', 'Edit Slide', 'Delete Slide' and 'Start Run' buttons become active



The Open Format i6000™ [Design Window] tool bar has the following options

Icon	Function
	<p>New File - Use this button to create a new staining run.</p> <p>Open File - Use this button to open a stored staining run.</p> <p>Save - Use this button to save changes to a staining run.</p> <p>The default filename is mm/dd/yyyy/rrr. Where 'm' is the current month (2 digits), 'd' is the current day (2 digits), 'y' is the current year (4 digits), and 'r' is the sequential run of the day (3 digits).</p> <p>After clicking Start Run, the following message appears:</p> <p>If you want to assign a unique filename to the staining run, click 'File' and select 'Save As', before clicking the 'Start Run' button.</p> <p>When the Stored Runs window appears, type the desired filename, and click Save or press the 'Enter' key. The filename can be any combination of letters and numbers.</p> <p>The file type, i.e. opn will be added automatically by the i6000™ software.</p>




Options- Use this button to display the Global Settings screen.

4.2 Global Settings

The user has the option to make various choices before the run starts. These choices will be used as the default settings for all slides within a run unless otherwise indicated.

Click on the '**Options**' button from the tool bar to open the '**Global Settings**' dialog box.

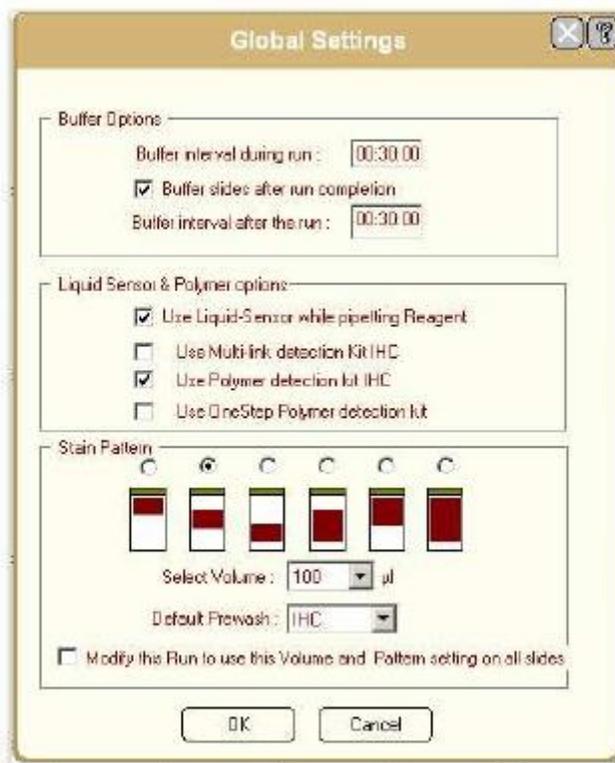

The following settings can be altered from this dialog box:

4.2.1 Buffer Options: The user can set a time to ensure that the instrument buffers the slide once within the chosen time interval during a run or after run completion.

4.2.2 Liquid Sensor Option: The user can run with or without the liquid level sensor.

4.2.3 Use Polymer Detection Kit for IHC: The user can choose polymer detection kit and protocol for IHC

4.2.4 Stain Pattern: The user can choose the staining pattern and default volume for all slides.



Check the box on the left side of the message “Modify this Run to use this volume & pattern settings on all slides”, to override all previous settings by the global settings.

	<p>Cut –Use this button to cut from the slide map one or multiple slides. Note that this button is only available if a slide has been selected.</p> <p>Copy – Use this button to copy contents of one slide to the clipboard. Note that this button is available only when one slide with protocol is highlighted.</p> <p>Paste – Use this button to paste from the clipboard to one or multiple slides. Note that this button is available only after a slide has been copied.</p> <p>Print Slides – Use this button to print the slide map</p> <p>Insert – Use this button to insert a protocol into the selected slide position(s). Note that this is only available when one or more slides positions have been selected.</p> <p>Edit – Use this button to edit one or multiple selected slides. Note that this button is only available when selected slides have a protocol inserted.</p> <p>Delete - Use this button to delete one or more slides from the Slide Map. Select the slides to be removed and click Delete. A confirmation message “Do you wish to Delete the selected Slide(s)?” will be displayed. Click Yes to delete, No to cancel and return to the Slide Map.</p> <p>Start - Use this button once you have completed the slide setup and would like to start the run.</p> <p>About – Use this button to display the About screen. This screen contains the software version of the i6000™ software, and current user information.</p>

5.0 Executing a Run

5.1 Inserting Protocols


- Place the slides in the slide racks and select the corresponding positions in the slide map of the Design Window. The selected slides will be highlighted.
- Right click on the target slide and select 'Insert Protocol' from the drop down list. Alternatively double click on the slide position to open the 'Insert Protocol' window.

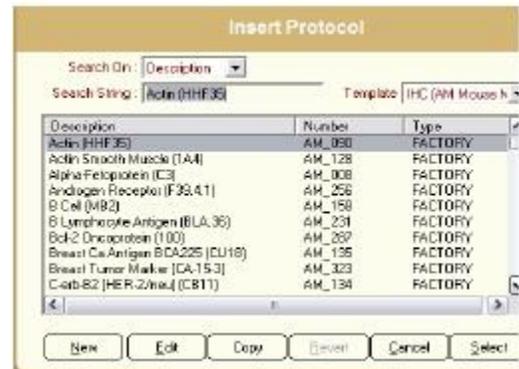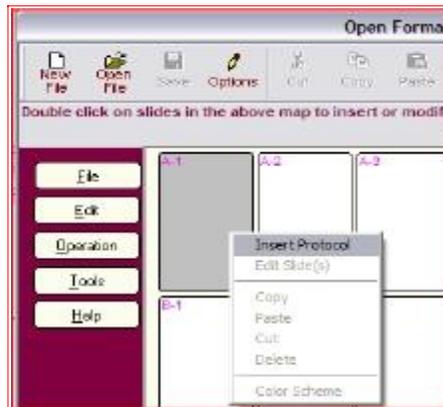

- The available protocols are displayed in alphabetical order along with the Catalog Number and Protocol Type.


Note: The antibody protocols currently listed in the software are provided as a reference for use with BioGenex primary antibodies and detection reagents only, and are subject to change. Please consult the appropriate BioGenex product inserts for the most current information before setting up a staining run.

- The user can search for the desired protocol using Description, Number, or Type as search criteria.

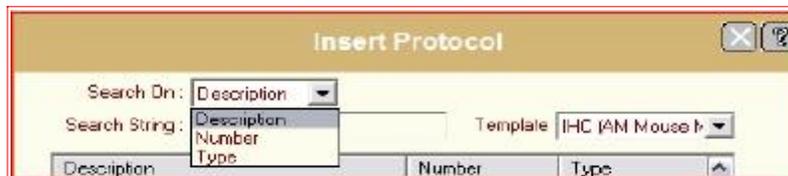

- Click the 'Search On' field to change the search criteria. Select the desired criteria from the drop down list. The same information should be entered in the 'Search String' field. The desired protocol is highlighted.
- Alternatively, use the scroll bar on to the right, to quickly select the protocol of interest from the list.



- Once the protocol is selected, click the 'Select' button to insert it into the target slide.

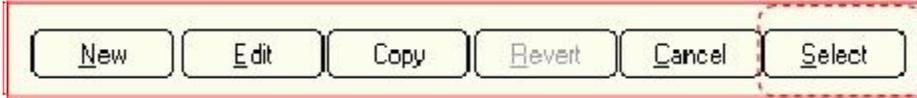

- Click the 'Cancel' button will return you to the i6000™ System [Design Window] without inserting the selected protocol.
- After inserting at least one protocol in any slide position, the 'Start Run' button will be enabled.


5.2 Copying Protocols

- To add the same protocol to multiple slides click 'Copy' on the tool bar, or right click on the slide with the desired protocol and select 'Copy' from the drop down list.
- Select the slide position(s) that will receive the desired protocol. To select consecutive slide positions, press the 'Shift' key while selecting. In case of random slide positions, use the Control key.
- Select 'Paste' from the tool bar, or right click on the target slides and select 'Paste' from the drop down list to insert the copied protocol.

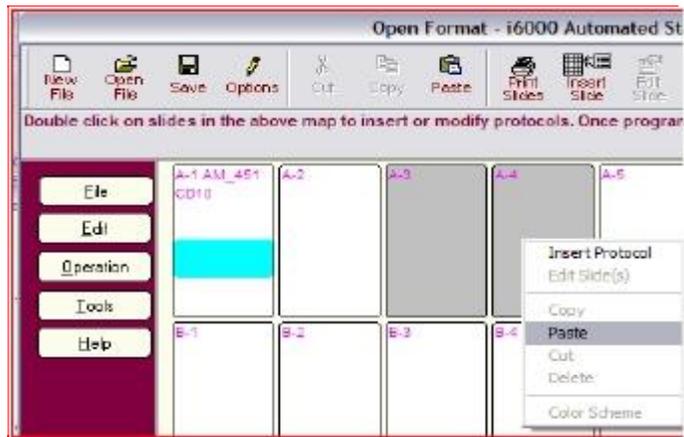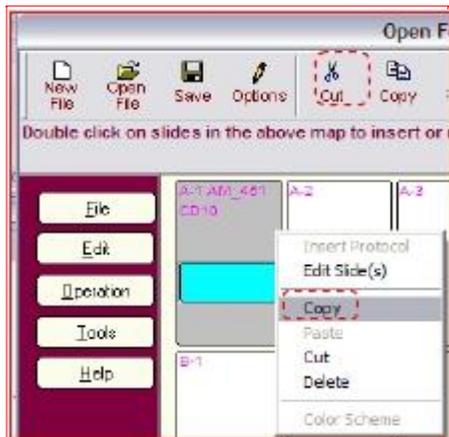

5.3 Edit Slide Parameters

Once the protocols are inserted into the slides, the user can edit the slide parameters as required. In order to edit, click once on the slide for its selection, and then right click on the same position to choose 'Edit Slide' from the drop down list. The **Edit Slide** window opens (Refer to section 6.3).

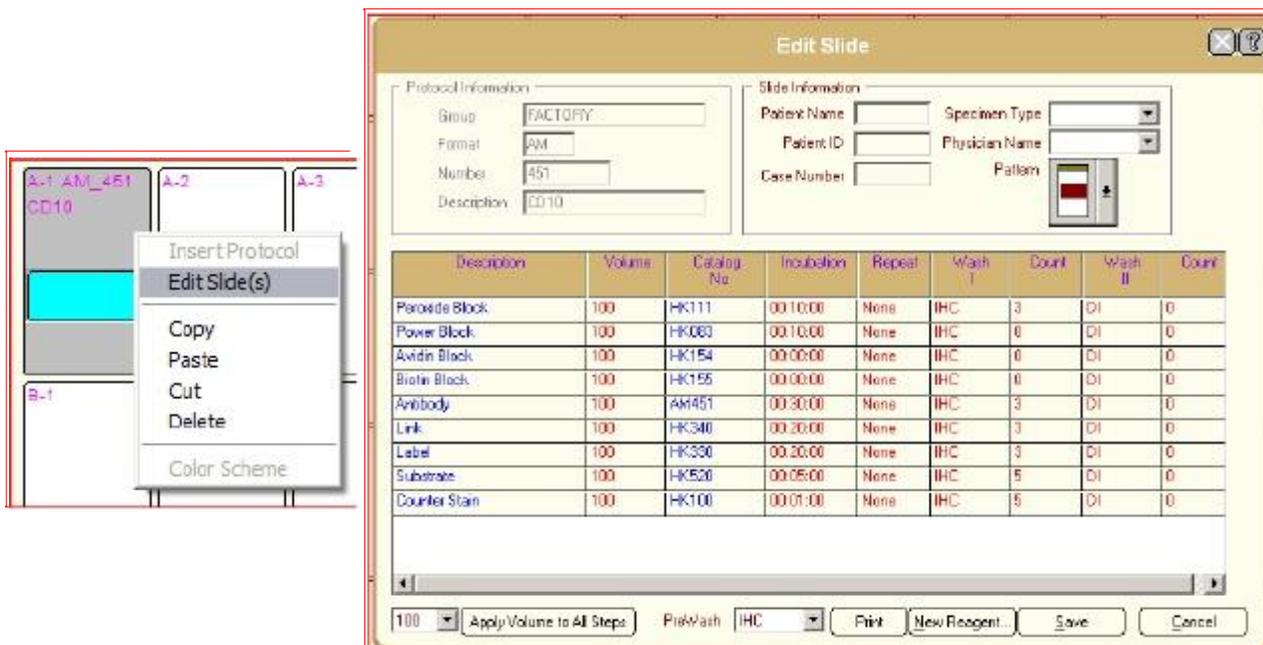

The Edit Slide screen is composed of data from the selected protocol and includes slide specific information. Slide specific information is listed on top right-hand side box titled Slide Information.

Note: Changes made in these steps will only affect the next run. The changes have no effect on the stored protocol.

- In order to make changes to the protocol, you must edit the protocol itself from the Edit Protocols window under 'Operation' Menu.
- Physician Names and Specimen Type are listed in the drop-down menu. Users may select from one of the listed names.
- If the desired name is not within the list, the user can add it.
- When finished editing parameters, click 'Save'. Selecting Save will return you to the Slide Map. Any editing of the displayed parameters will be saved only for that slide position.


5.4 Deleting Slides

Any slide or group of slides can be deleted.

- Select the target slide(s) to highlight and right click.
- Click on 'Delete' from the drop down menu.



- A confirmation dialog box will be displayed on the screen. Click 'Yes' to delete the slides.

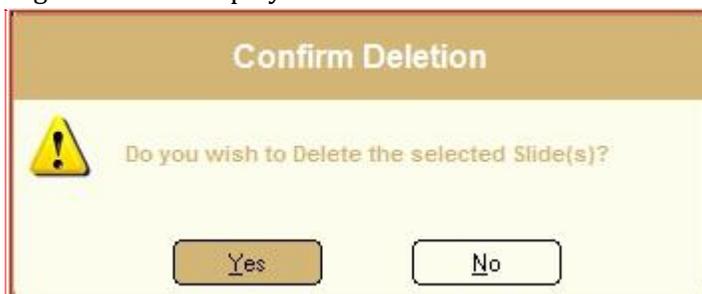

Choosing 'No' will take you back to the i6000™ System [Design Window] without deleting the selected slides. Slides can also be deleted by clicking on the Cut button.

5.5 Starting the Run

- Once the protocols are inserted, click the '**Start**' button to initiate the staining run.
- The run will be automatically saved with a default filename and file type before the run begins.

Note:

- The default filename is **mmddyyyyrrr**, where 'm' is the month, 'd' is the current day, 'y' is the current year, and 'r' is the run sequence number of the day.

- To assign a unique filename to the staining run, select **File> Save As**, before clicking the Start button.



- Click 'OK' to open another dialog box




- To insert STAT or continuous slides, and to save time and to shorten the time of run interruption, you could start your initial run with enough volume to accommodate subsequent slides. However, this is not necessary as you can always insert the required reagents at a later time.
- Press 'OK' to display the reagent screen.

5.6 Reagent Screen

This screen displays the positions in the reagent rack where you must load each reagent with sufficient volume. The catalog number, description, and minimum volume required for all reagents needed to complete the staining run are displayed on the Reagent Screen.

- To print the Reagent Map click on the Print button on the top right corner.
- Once all the required reagents are loaded, click on the Proceed button to continue the run. The 'Load Consumables' screen will be displayed.

- Ensure that all the consumables needed for the Run are loaded in the system. Verify and Check the box against each item to proceed with the run.
- Click 'Proceed' button to open the 'Prime Solutions' dialog box

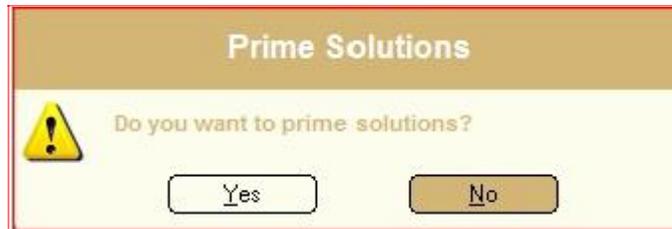

- Click 'Yes' to open the 'Prime Solution' window.



- Specify the time in the "Prime for – secs" field. Check the boxes against the solutions to be primed. Click the 'Start priming' button.
- Prime the solutions again if required. Once done, click the 'Close' button to open the 'specify the pipette position' screen.



- Select an appropriate option and click 'Continue'.

5.7 Delayed Run

- The system prompts the user whether to start a '**Delayed Run**' or not? Click '**No**' to start the run right away.



- To start a delayed run, click 'Yes'. The following screen will be displayed



- The current date and time are the default entries. Enter the date and time that you want the run to start (maximum delay is 3 days). You can access a calendar by clicking on the arrow to the right of the date.
Note that the time is entered in 24-hour format.
- You can buffer the slides while waiting for the run to begin by selecting the Buffer Slides during waiting option. The frequency of buffering is set in the global settings (refer to the global setting and other tools section).
- Clicking Cancel will start the run with no delay. Click Start Delayed Run to accept the time entered,
- The 'Delayed Run Timer' window will be displayed.




- You must leave the system on and the program running. Closing the program will cancel the programmed delayed start run even if the program is restarted before the selected start date and time.
- There is a countdown window that appears until the run begins.
- This window gives you the option to cancel the run by clicking on the Stop Run button, or start the run immediately by clicking on the 'Start Now' button.

- The number of buffering cycles completed on each slide during the waiting period is also displayed.

5.8 Slide Information

Once a run is started, the following information will be displayed for each slide:



- Slide position
- Antibody number
- Antibody description
- Protocol Step that is currently being performed
- Remaining Time for that Step (countdown)
- Staining Pattern (colored area)
- Current step number and total steps of the protocol

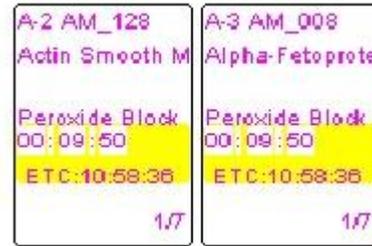

5.9 Slide Priority

There are two priority settings for any slide position, 'Normal' and 'High'. Initially all the slides are at normal priority by default. The user can prioritize the slides as required in order to process them on a higher priority than the remaining slides in that run .

- To increase the priority of a slide, right click on the desired slide position and from the drop down list, select '**Increase Priority**'. This option is available only for normal priority slide(s)
- To decrease the priority, right click on the same position and select '**Decrease Priority**' from the drop down list. This option is available only for high priority slide(s).

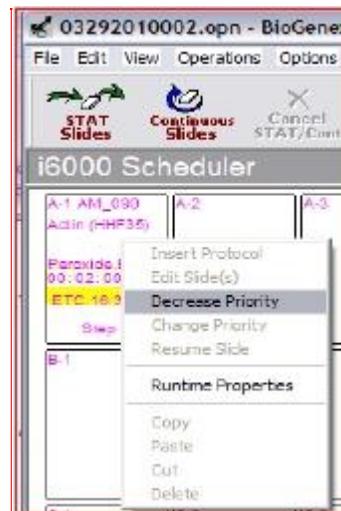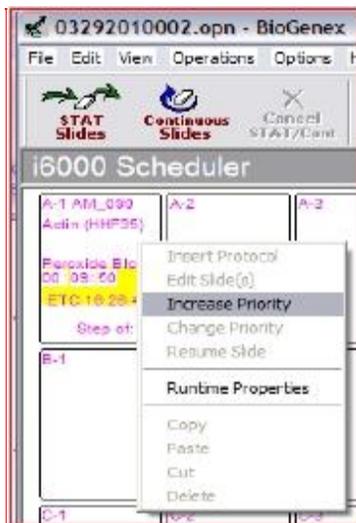

5.10 Run Information

The user can view the estimated time for completion of an ongoing run. The run completion time can be influenced by a number of factors like buffering cycle, addition of STAT or continuous slides, fixing errors like dropped pipette tip or lack of reagents. The system recalculates the run completion time after each of these steps.


To view the run completion time, select the **'View > Run Information'** to open the **'Run time estimation'**

Dialog box with a message "Run is expected to complete at ---".



5.11 Runtime Properties


The status of each slide can be easily tracked from the information displayed on the respective positions. More detailed information can be obtained from the **'Runtime Properties'**. To find out individual slide progress, right click on the slide and select **'Runtime Properties'** from the drop down list to open the **'Properties of Slide [-]'** window. It displays the information on the steps that are ongoing, completed or yet to begin.

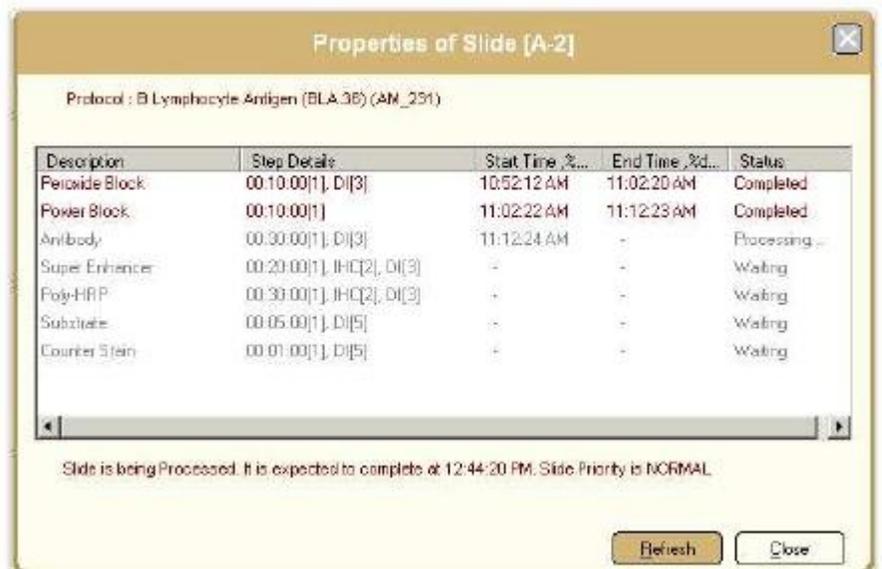

5.12 Inserting Continuous Slides

Initiating a Continuous Run - The i6000™ system facilitates optimal workflow through the continuous access. Depending on the availability of free racks, the user can add extra slides to an ongoing run, increase priority and save time.


- To start a continuous run, click on the '**Continuous slides**' button.
- The ongoing run is paused for 2min with a message "**CONTINUOUS enabled. Insert slides and click on 'Continue Run' button...**" and a '**2 min Timer**' displayed below the tool bar

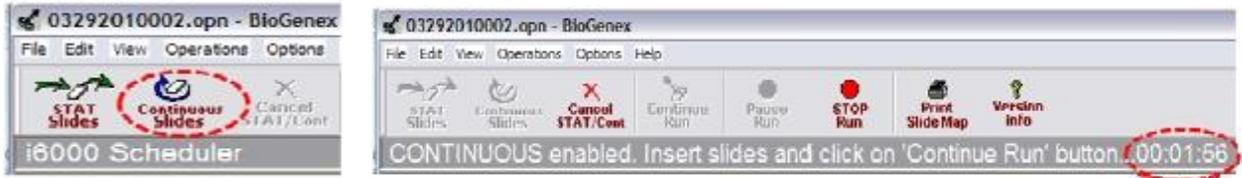

- Insert the new slide(s) into the vacant positions available.
- Insert the protocol by double clicking on the target positions to open the 'Insert Protocol' window.
- Select the protocol of interest from the list available, and click the 'Select' button to insert the protocol into the target slide position.
- Click on 'Continue Run' to open the following message.



- Click '**OK**' to continue with the run.
- If additional reagents are required, the 'Reagent Screen' will appear.
- Load the reagents in the same position specified in the reagent screen and click on 'Proceed' to continue the run.
- A delay exceeding 2min. in initiating a continuous run sets off a system alarm along with a warning message <**Timeout Occurred**> displayed on the screen.
- Click 'No' to reset the timer and continue setting up the continuous run. Click 'Yes' to cancel the continuous run.

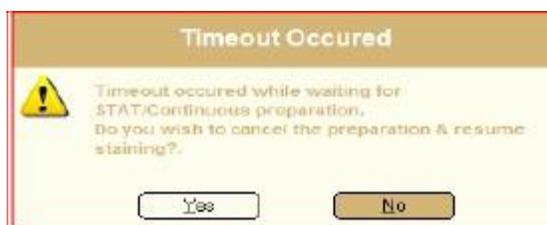

It is advisable to have all the reagents and consumables ready before initiating a continuous run.

Note: The user can attempt only thrice to add the continuous slides and their number should not exceed 30. If the limit is exceeded, the following dialog box will be displayed on the screen. Click 'OK' to continue with the ongoing run.




5.13 Inserting STAT Slides

STAT access is a very handy feature that allows the user to add a maximum of six slides to an ongoing run (depending on the availability), which are processed at a higher priority than normal.

- Click on the <STAT Slides> button in the tool bar to start the STAT run.



- The ongoing run is paused for 2min with a message “**STAT enabled. Insert slides and click on 'Continue Run' button...**” and a **2 min Timer** displayed below the tool bar.
- Place the new slide(s) into the vacant positions available.
- Insert the protocol by double clicking on the target positions to open the 'Insert Protocol' window.
- Select the protocol of interest from the list available, and click the 'Select' button to insert the protocol into the target slide position.
- Click on 'Continue Run' to proceed ahead. If the reagents for the new slide(s) are already available with the system, a message stating the same will be displayed.

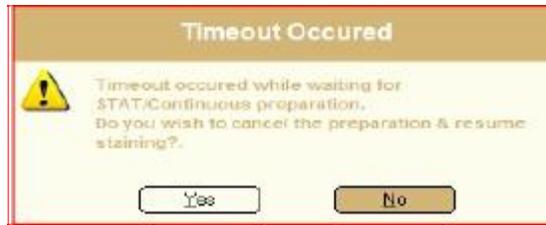

It is advisable to have all the reagents and consumables ready before initiating a continuous run.

5.14 Run Completion

At the end of the run, the **'Run Completed successfully'** dialog box is displayed along with the following information

- I. The time of run completion
- II. How long ago the run completed
- III. The number of buffering cycles applied since the run ended

- Click on 'Print Reports' button to open the 'Welcome to Reports i6000' dialog box.
- On the left panel of this dialog box, following options are available:
 - Run Summary Report
 - Case Summary Report
 - Run Log Summary
 - Assay Report
 - Reagent Summary by Date
 - Reagent Summary by catalog
 - Reagent Detail Report



Once an option is selected, a brief description of its function is displayed on the right panel. Select an option and click **'Next'** button to open the report in a new window in a preset format.


6.0 Protocols

A protocol is a set of operations performed sequentially under defined conditions. The i6000™ Infinity system comes along with 'FACTORY' protocols that are optimized to deliver consistent and reliable results when using BioGenex reagents. However, the user can modify the factory protocol as per the requirements or create a new protocol. If a **FACTORY** protocol is modified, then it is saved as **FACTORY_USER** protocol. A new protocol created by the user will be saved as 'USER' protocol.


6.1 Creating a New Protocol

All the protocols created by the user are considered USER protocols. To create a new protocol:


- Click on Operation from the i6000™ System [Design Window] menu, and select Edit Protocols from the drop-down list.
- The 'Edit Protocols' screen will be displayed.

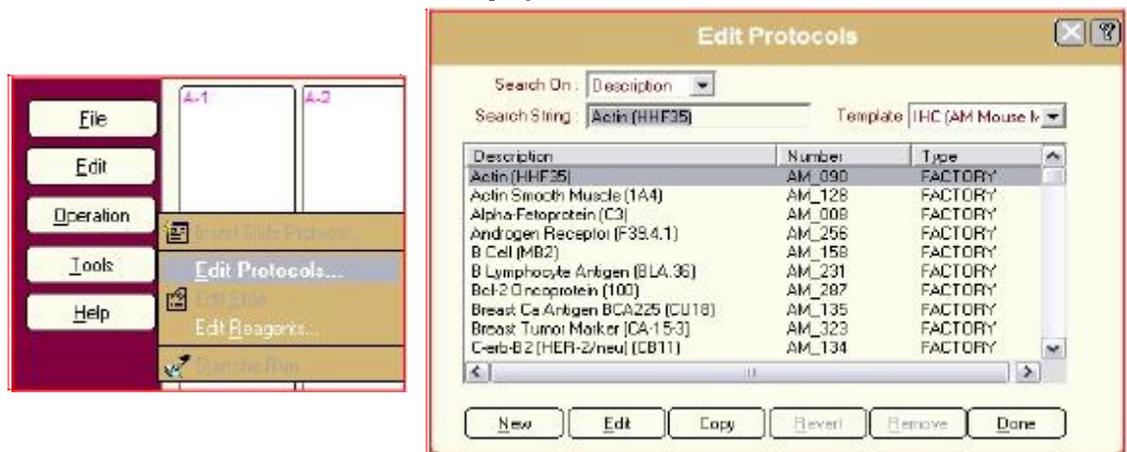

- Click on the 'New' button to open the 'Research Steps' dialog box. Specify the number of steps in the box provided and click OK.
- The 'Create a New Protocol' window will be displayed.

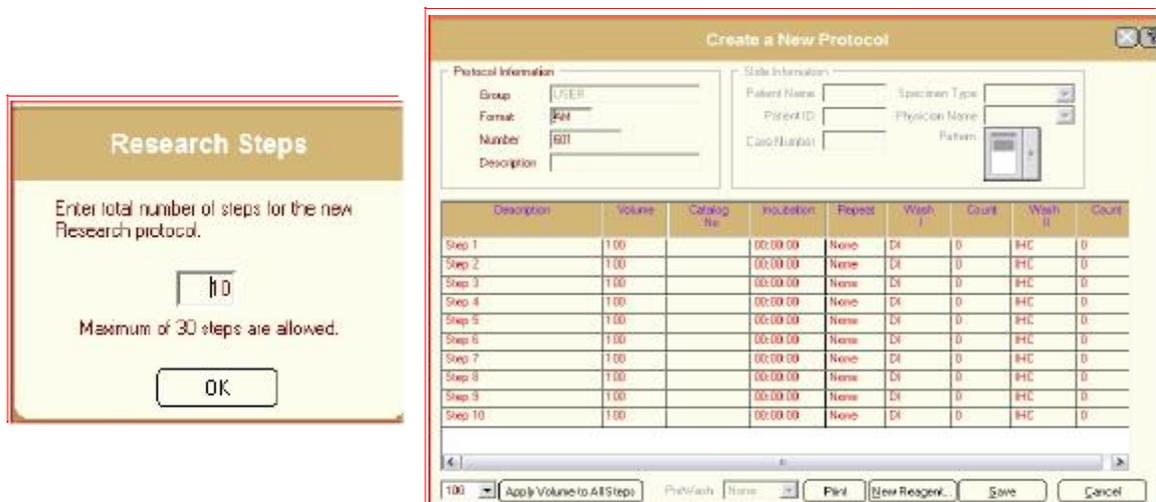


- In the Protocol Information area, type in the antibody number of your choice starting from 601 to 999, and your antibody description.
- All the steps (highlighted in red) can be edited.
- Each step can be renamed as follows:
 - Click on the corresponding cell to visualize an 'Arrow' button by the side of it.
 - Click the arrow to select the appropriate name from the drop down list.
 - Click on the name to assign it to the step

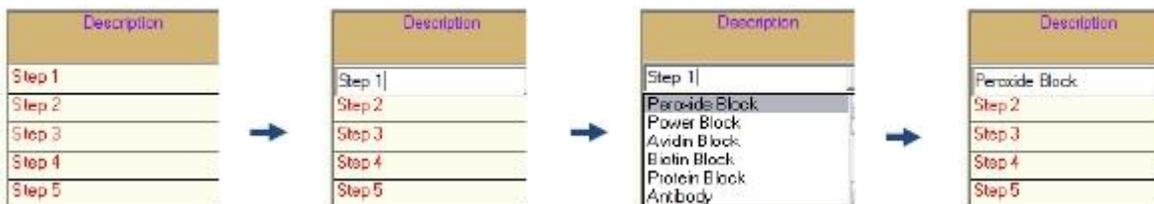


- Similarly to edit a parameter, click on the corresponding cell to visualize an 'Arrow' button by the side of it.
- Click the arrow to open a drop down list.
- Click on the appropriate value to assign it to the cell.

If the list of Reagents provided does not contain what you want to use, add a new reagent to the database by clicking on the 'New Reagent' button. The following dialogue box will be displayed.

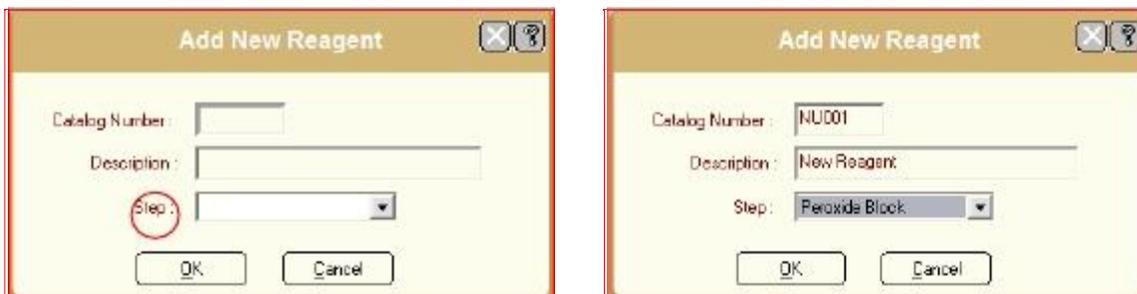

- Fill in the information in the respective fields. The format for entering the Catalog No. is 'XX###', where 'X' is any letter and '#' any digit between 0-9.
- To name a step, click on the arrow button and from the drop down list, select the appropriate name.
- Click 'OK' to add the new reagent to the database. The 'Create a New Protocol' window is displayed on the screen.
- To select the newly added reagent, click on the corresponding step in the 'Catalog No.' column and select it from the drop down list.

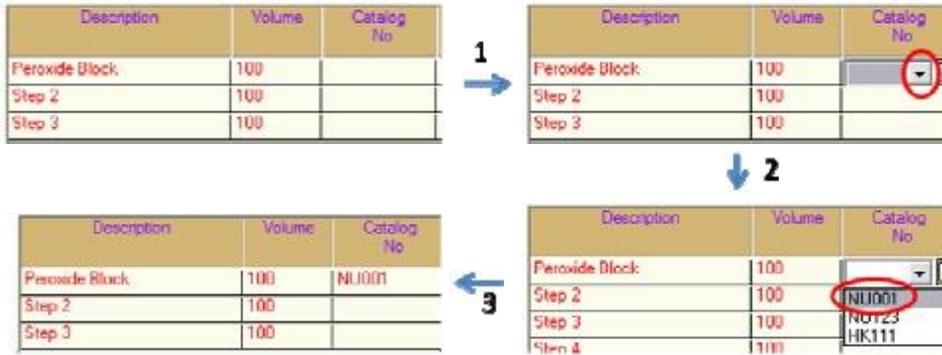


- When all the steps are incorporated, click on 'Save' to return to the 'Edit Protocols' window. The new protocol is saved as USER type and is highlighted.
- Clicking 'Cancel' will return you to the 'Edit Protocols' window without creating the USER protocol.


Creating a New Protocol Using an Existing Factory Protocol as a Template

6.2 Copying a Protocol

- Click on Operation from the i6000™ System [Design Window] menu, and select Edit Protocols from the drop-down list.
- The 'Edit Protocols' screen will be displayed.

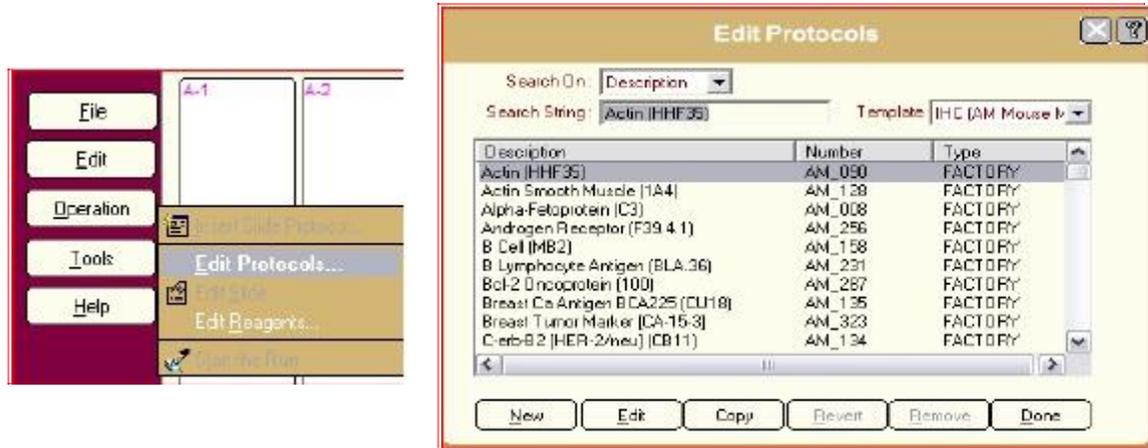

- Select a 'FACTORY' protocol of interest and click the 'Copy' button to open the 'Create a New Protocol' window.

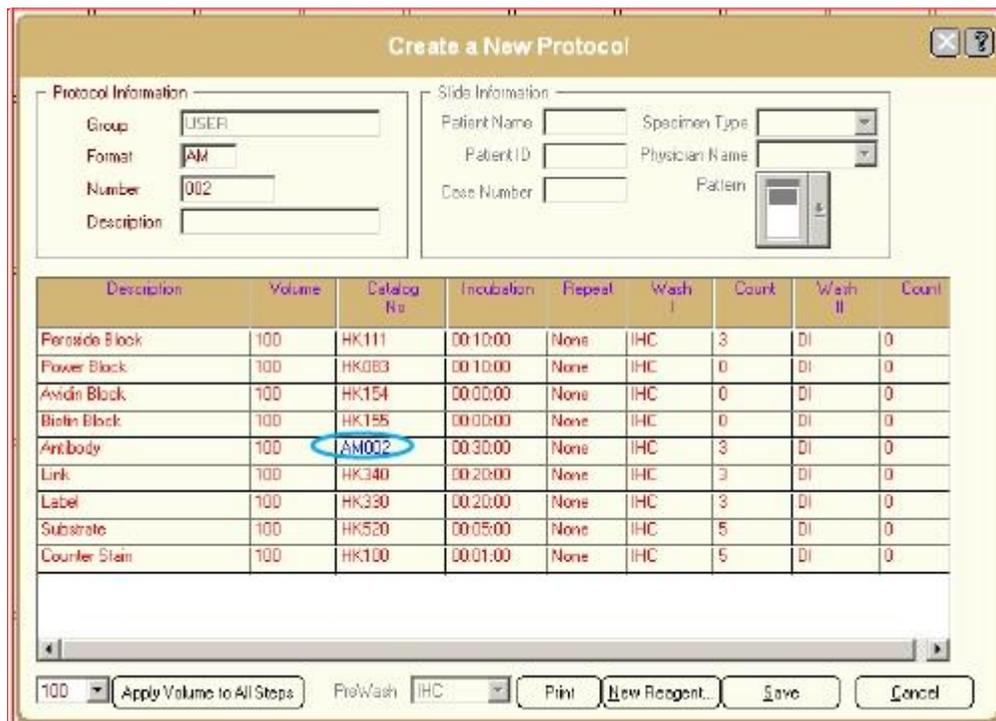


- All the steps can be modified as described in the previous section, except for the antibody (highlighted in blue).
- When all the necessary modifications are done, click 'Save' to return to the 'Edit Protocols' window. The new protocol is saved as USER type and is highlighted.

- Clicking 'Cancel' will return you to the 'Edit Protocols' window without creating the USER protocol.

6.3 Edit Protocol

A 'FACTORY' protocol when edited, changes to a '**FACTORY_USER**' protocol. Only the Administrator is authorized to edit a FACTORY protocol. To edit a Factory protocol:

- Open the <i6000 Automated Staining System [Design Window]>
- Click on the '**Operation**' menu and select '**Edit Protocols**' from the drop down list.

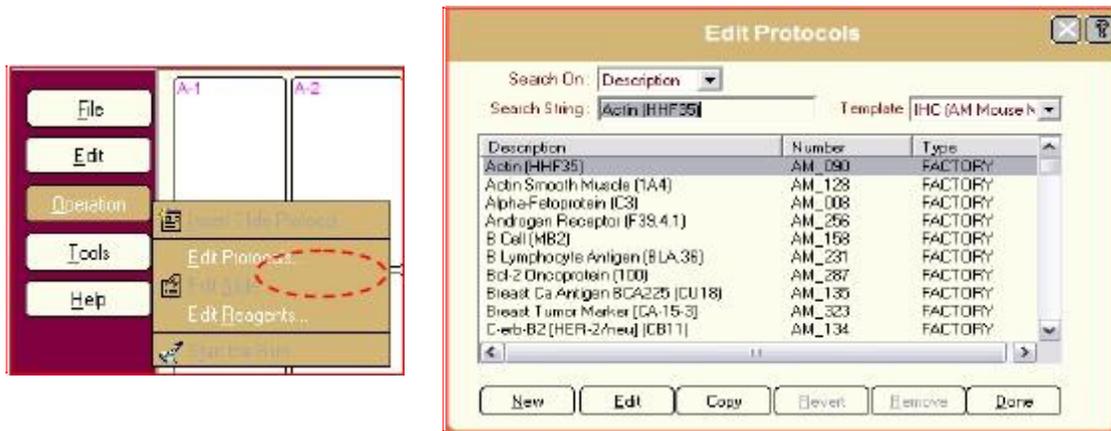

- The '**Edit Protocols**' window will open. Select the protocol of interest from the list available, by clicking on it.
- Click the '**Edit**' button, or double click on the highlighted protocol.
- A warning message "Editing a FACTORY protocol will result in a FACTORY_USER protocol. Do you wish to continue?" will appear on the screen.
-



- Click '**Yes**' or press '**Enter**' key to open the '**Edit Protocol**' window. Click '**No**' to return to '**Edit Protocols**' window

Edit Protocol

Protocol Information

Group:

Format:

Number:

Description:

Slide Information

Patient Name:

Patient ID:

Case Number:

Specimen Type:


Physician Name:

Patient: 

Description	Volume	Catalog No	Incubation	Repeat	Wash I	Count	Wash II	Count
Peroxide Block	100	HK111	00:10:00	None	IHC	3	DI	0
Power Block	100	HK083	00:10:00	None	IHC	0	DI	0
Avidin Block	100	HK154	00:00:00	None	IHC	0	DI	0
Biotin Block	100	HK155	00:00:00	None	IHC	0	DI	0
Antibody	100	AM090	00:30:00	None	IHC	3	DI	0
Link	100	HK340	00:20:00	None	IHC	3	DI	0
Label	100	HK330	00:20:00	None	IHC	3	DI	0
Substrate	100	HK520	00:05:00	None	IHC	5	DI	0
Counter Stain	100	HK100	00:01:00	None	IHC	5	DI	0

PreWash:

- All the parameters can be edited except for those items in blue.
- To edit a parameter, click on the corresponding cell to visualize an 'Arrow' button by the side of it.
- Click the arrow to select from the drop down list. Click on the appropriate value to assign it to the cell.

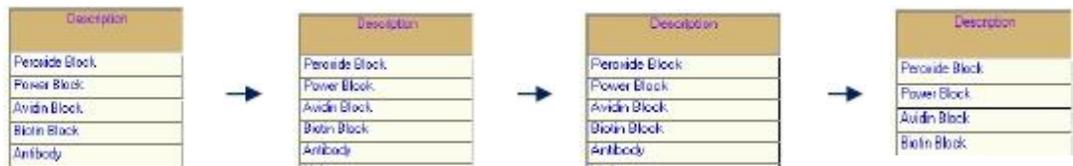


- Click 'Save' button to save the changes as **FACTORY_USER** Protocol. This will be highlighted when you return to the 'Edit Protocols' screen.
- Click 'Cancel' button to undo all the modifications and return to the 'Edit Protocols' screen.

6.4 Revert Protocol

A **FACTORY_USER** Protocol can always be reverted back to the original **FACTORY** protocol settings from which it is derived. Only the Administrator is authorized to revert back a **FACTORY_USER** protocol to a **FACTORY** protocol. To revert a **FACTORY_USER** protocol

- From the Operation menu click on 'Edit Protocol' to open the Edit Protocols window.
- Select the **FACTORY_USER** protocol that needs to be reverted. Once selected, the Revert button is activated.
- Click the Revert button. A 'Revert Operation' dialog box is displayed on the screen.

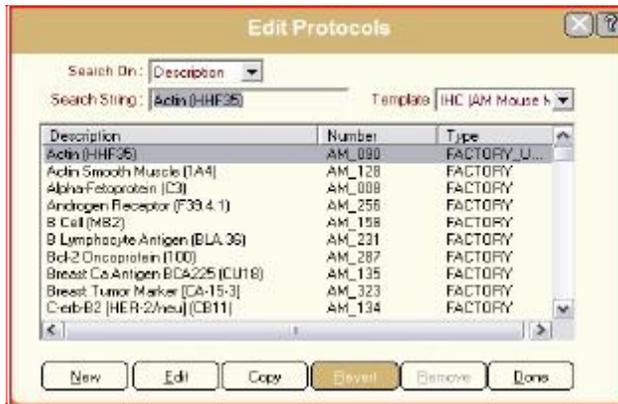

- Click 'Yes' to restore the Factory Protocol settings. Once restored the 'Edit Protocols' window will be displayed again.
- Click 'No' to retain the FACTORY_USER protocol.

6.5 Deleting a Protocol

- Only USER protocols can be deleted. To delete a User protocol,
- Open the 'Edit Protocols' window from the 'Operation' menu.
- Select the USER protocol to be deleted.
- Once any USER protocol is selected, the 'Remove' button gets activated.
- Click on the 'Remove' button. A message "This operation will delete the protocol 'XX_###'. Do you wish to continue?"
- Click 'Yes' or press the Enter key to delete the protocol.



Button Functions in 'Edit Protocols' window

New:	Click this button to create a new protocol
Edit:	Click this button to open the 'Edit Protocol' window
Copy:	Click this button to copy all the steps involved in a protocol of interest. Once copied, the steps can be modified and saved as USER protocol

Revert:	Click this button to undo all the changes in a FACTORY_USER protocol and restore the original FACTORY protocol. This button is active only when a FACTORY_USER protocol is selected
Remove :	Click this button to delete a USER protocol from the data base. This button is active only when a USER protocol is selected.
Done:	Click this button to close the 'Edit Protocol' window

Button Functions in '**Edit Protocol**' window


Apply Volume to All Steps:	Click this button to apply the volume specified in the adjacent box to all the steps in the protocol. To set the volume, click the arrow button and select an appropriate value from the drop down list.
PreWash	Click this button to select IHC or DI buffers to prewash the slides.
Print:	Click this button to print the protocol.
New Reagent:	<p>Click this button to open the 'Add New Reagent' dialog box. There are three fields in this dialog box, viz., Catalog Number, Description, and Step.</p> <ul style="list-style-type: none"> • Type in the required information as applicable in the Catalog Number and Description fields. • Click the arrow button in the 'Step' field and select the appropriate step from the dropdown list. • Click 'OK' to add the new reagent to the database. Click 'Cancel' to close the dialog box without adding the new reagent to the database
Save:	Click this button to save the changes made to the protocol.
Cancel:	Click this button to close the window without saving the changes.

7.0 Reagents

7.1 Creating a New Reagent

All reagents created by the user are considered USER reagents.

- To create user-defined reagents, click on 'Operation' from the i6000™ System [Design Window], and select 'Edit reagents' from the drop-down list.
- The 'Reagent Database' window is displayed on the screen.

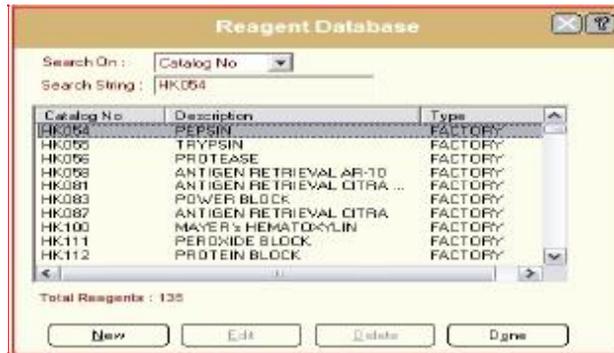

- Click 'New' to open the 'Add New Reagents' window. There are three fields in this window, viz., Catalog Number, Description, and Step.



- Fill in the information in the respective fields. The format for entering the Catalog No. is 'XX###', where 'X' is any letter and '#' any digit between 0-9.
- To name a step, click on the arrow button and from the drop down list, select the appropriate name.

7.2 Editing a Reagent Description

Only USER reagents can be edited. Click on Edit reagents from the Operation drop down menu. The following screen will be displayed.

8.0 Global Settings and Other Tools

8.1 User Manager

When the software is first installed, it can be accessed only by the Administrator. The 'User Profiles' need to be added by the administrator. To access the User Manager:

- Select the **Tools** menu and click on **User Manager** from the drop down list. The User Manager screen will be displayed.



- Click the 'New' button to display the 'New User' window
- Enter a User ID, the full Name, Password, and Confirm Password. The entry made in 'Confirm Password' should match exactly with that of the 'Password'.
- Press 'OK' to save the new user. The user will be added to the user database.




8.2 Edit User Profiles

- Select a user from the User Manager window and click the Edit button.
- The **'Edit User'** window will be displayed.
- Make the necessary changes and click **'OK'** to save.



8.3 Delete User Profile

- Open the User Manager window and select the user profile to be deleted. The 'Remove' button gets activated only when a user profile other than the Administrator is selected
- Click '**Remove**' button to open a confirmation message.

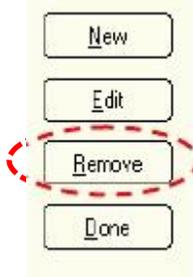

- Click 'Yes' to delete the user profile from the database. The user will no longer have any access to the system.
Note: Administrator profile cannot be deleted.

8.4 Change Password

While creating a user profile, the Administrator will assign a password. The user can log in using this password. However, it can be changed if the user intends to do so.


- From the 'Tools' menu, click on 'Change Password'.
- The 'Change Password' window will be displayed.
- Enter the new password and retype the same in 'Confirm Password' field.
- Click 'OK' to save the new password. Clicking cancel will close the window without saving the changes.



8.5 Physician Names

A list of Physicians Names can be maintained for easy entry in the **Edit Slide** window. Each name entered will be available in a drop down list in the Slide Information area of 'Edit Slide' Window. To enter Physician names:

- From the **'Tools'** menu, click on the **'Physician Names'** to open the **'Physician Entry Form'**
- Type the name in the **Name** field and click **Add**. The name will appear on the list below.
- To delete a name from the list, select the name and click **Delete**.

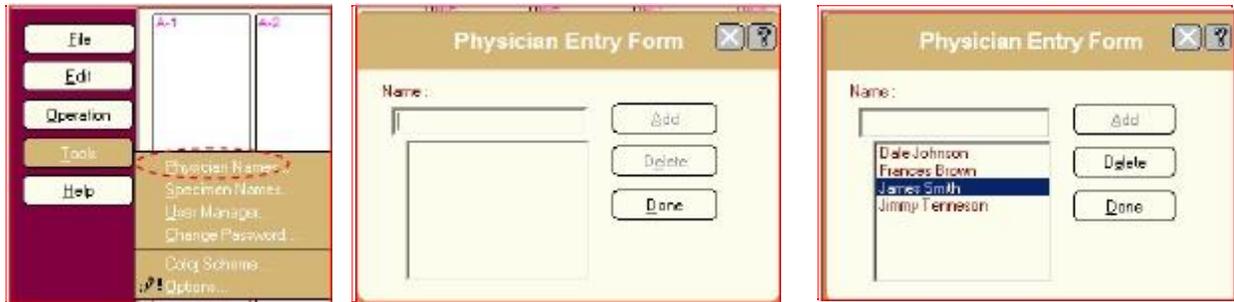

8.6 Specimen Types

A list of specimen types can be maintained for easy entry in the 'Edit Slide' window. Each name entered will be available in the drop down list in the Slide Information area of 'Edit Slide' Window. To add the specimen names:

- From the **'Tools'** menu, click on the **Specimen Names** to open the **'Specimen Entry Form'**
- Type the name in the **Name** field and click **Add**. The name will appear on the list below.
- To delete a name from the list, select the name and click **Delete**.



8.7 Color Scheme

The color scheme is meant for representing the status of the following:

- Slides not started
- Slides in progress
- Completed slides

- Slides with Error
- Text on the slide

From the **'Tools'** menu, select **'Color Scheme'** to open the **'Color Scheme'** window



9.0 Pause or Stop Function

Prolonged pause during a staining run may cause invalid test results due to extended incubation time for the affected step and the slide(s).

9.1 Emergency Stop

The automated system has an Emergency off Switch, located on the front side of the instrument in the lower-right hand corner. When pressed, the 'System Error' screen is displayed.





To continue, close the application by clicking on the **Close** button, turn back the Emergency stop and start a new run.

9.2 Software Pause

To pause a run at any time

- Click on the **'Operations'** button in the menu bar and select **'Pause'** from the drop down list.

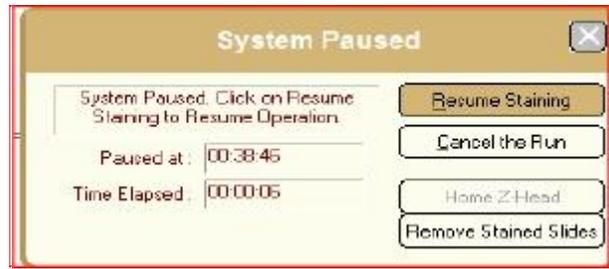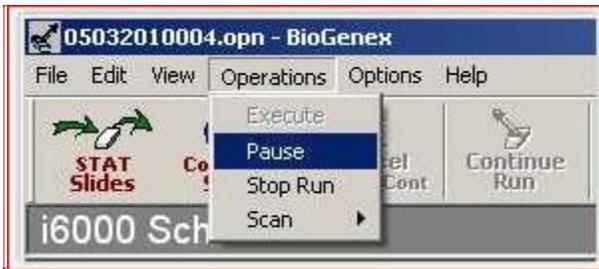

- Click the '**Resume Staining**' button or press '**Enter**' key to continue the staining run.
- Click '**Cancel the Run**' to cancel the run
- Click '**Home Z-Head**' to return the Z-head to the home position
- Click '**Remove Stained Slides**' to remove the completed slides from the slide map.

9.3 Software Stop

- To stop a current run, click on the red **Stop button**, or select 'Stop Run' from the **Operations** menu.
- The '**Stop Confirmation Message**' will be displayed



- Click '**No**' to return to the staining run.
- Click '**Yes**' to terminate the current run. The **Reset Machine** dialog box will be displayed
- To keep the instrument pressurized, click on the 'No' button. To release the system pressure, click on the 'Yes' button.



10.0 Reports

The i6000™ software provides the user with an array of reports summarizing data for easy reference. Users can design and customize new reports.

10.1 Starting the Reports Module:

There are two ways to launch the Reports module –

- 1) Print Reports from main log in screen- **“Print Reports”** window
 - 2) Once the staining run is completed- **“Run Completed Successfully”** window
- Click **‘Print Reports’** from the main log in ‘Welcome to i6000’ screen or the ‘Run Completed Successfully’ window



- The **‘Welcome to Reports i6000’** window is displayed



- Select a report format from the list available by clicking on it. A brief description of the report format selected will be displayed on the panel onto the right.
- Select a report format of interest and click **‘Next’** to proceed ahead with the print.
- Click **‘Cancel’** to go back to the main log in window.

Report	Function
Run Summary Report	This report prints out details of all the runs for a chosen date or range of dates
Case Summary Report	This report prints out all the slides that were composed for a specific case number
Run Log Summary	This report prints out the details for a specific run
Assay Report	This report prints out protocol details of all slides for a specific run
Workload Report	This report lists out all the reagents used in a specific run
Reagent Summary by Date	This report prints all the reagents used in a range of dates. Report is grouped by date.
Reagent Summary by Catalog	This report prints all the reagents used in a range of dates. Report is grouped by Catalog number
Reagent Detail Report	This report prints out the detailed usage of a specific catalog

- After selecting a report format, clicking on the 'Next' button would open a dialog box that can be used to set the parameters like date, run number, assay etc.
- For each report the desired parameters can be selected from the drop down list displayed by clicking the arrow button adjacent to the fields.

- Once the parameters are set, click on '**Preview & Print Report**' button.

- The chosen report will be displayed on the screen. Click on the printer icon



Assay Report



BioGenex
ACCELERATING THE RACE OF PRECISION MEDICINE

Summary of Protocols for the Run Number : 06152010001.opn

Print Date/Time: 6/15/2010 - 19:37:36

Instrument ID : AS60XXX
User ID : Administrator

1

Protocol : AM_451 - CD10	Start Time : 19:29:50	Patient ID :
Status : Incomplete	End Time :	Physician Name :
Run Date : 6/15/2010		Specimen Type :
		Case Number :


Step	Description	Reagent	Incubation	Volume	Wash I	Wash II
1	Peroxide Block	HK111	00:10:00 [1]	100 µl	DHC [3]	-
2	Power Block	HK083	00:10:00 [1]	100 µl	-	-
5	Antibody	AM451	00:30:00 [1]	100 µl	DHC [3]	-
6	Link	HK340	00:20:00 [1]	100 µl	DHC [3]	-
7	Label	HK330	00:20:00 [1]	100 µl	DHC [3]	-
8	Substrate	HK320	00:05:00 [1]	100 µl	DHC [5]	-
9	Complete Stain	HK100	00:01:00 [1]	100 µl	DHC [5]	-

Staining Quality: _____


Comments: _____

10.2 Export a Report

- To export a report, select a specific report and click on the 'Export Report' button

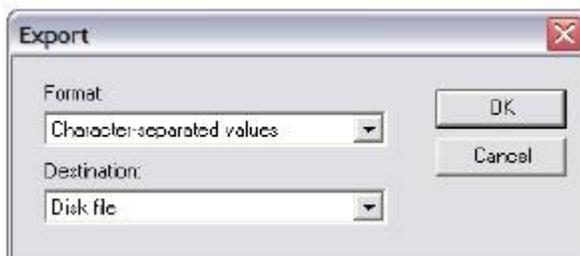

- Click the arrow buttons adjacent to the fields, and select the format and destination from the drop down list. The report file will be exported.

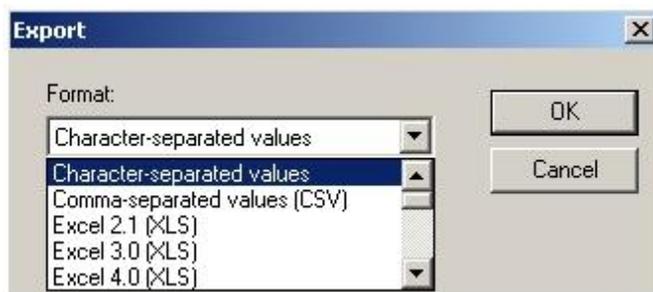


- There are different format options available for exporting a report. Click on the arrow adjacent to the Format option and select the format in which the report to be exported.

11.0 Calibration and Diagnostics

To start the calibration program, click 'Calibrate **i6000**' button from the 'Welcome to i6000' window. The 'Please make sure...' dialog box will be opened.



Press 'OK' to display the 'i6000 Calibration & Diagnostics Tool' window.

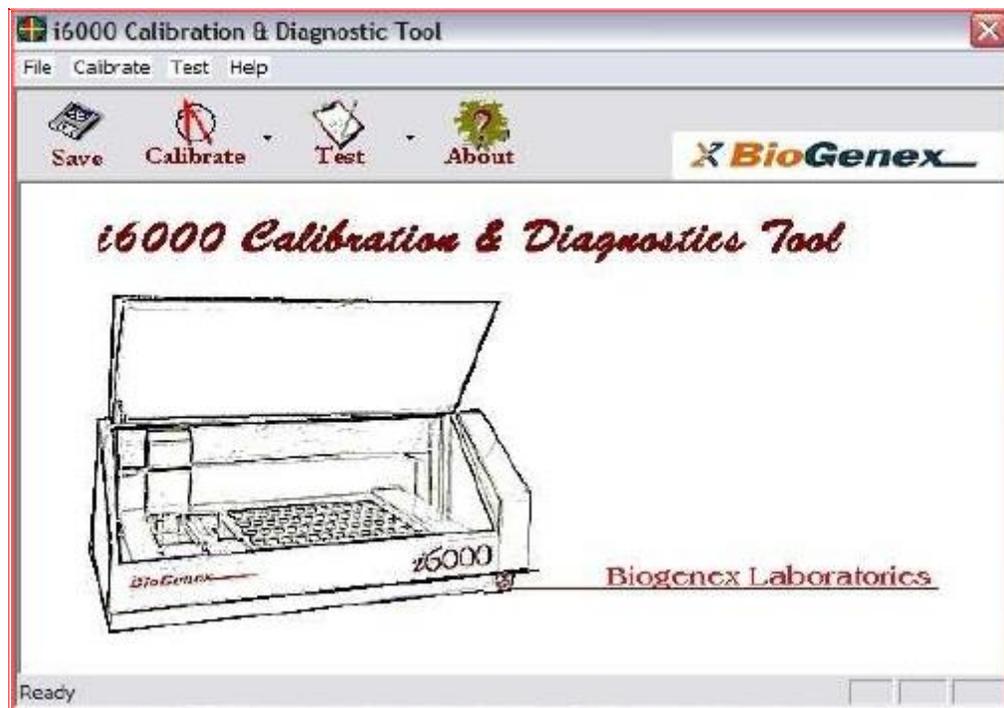


11.1 One Point Calibration

The i6000™ Infinity '1 Point Calibration' is made possible by the use of a one piece-working base. This base provides consistent spacing between functional locations of the system, i.e., the pipette tips, reagent racks, and slide racks are in fixed locations




- Select '1 Point Calibration' from the Calibrate menu.
- Insert the calibration pointer into the pipette tip adapter.
- The flat end of the pointer goes into the orifice at the bottom of the adapter. The pointer must be inserted completely; make sure that the flat disk portion of the pointer is touching the tip adapter.

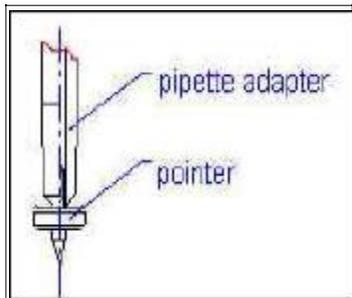

Caution: Do not move the Z-head when installing the pointer as this will cause inaccurate calibration.

- Once the pointer is installed, click 'OK' to display the 'Calibrate Positions' window.

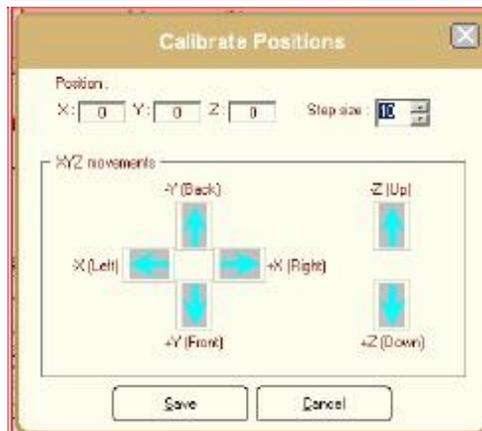


- Set the **Step Size** to 10.
 - Use the arrow buttons or the keyboard to move the Z-head to a position just above the calibration target.
 - Set the step size to 3
 - Lower the Z-head until the pointer is close to the center point of the calibration target. Adjust the X and Y positions again if needed.
- Caution:** Do not overdrive the pointer into the target. If you force the motor to 'skip', the calibration will be incorrect. Click '**Cancel**' and start the calibration all over again.
- Click on Save to move the Z-head from the target. Remove the calibration pointer to complete the One Point Calibration.

11.2 Multi-point Calibration

The tasks involved in this must be completed in the order they are presented in the Multi-point Calibration menu.

- Select Calibrate from the i6000™ Calibration & Diagnostic Tool menu. Choose '**Multi Point Calibration**' from the drop down menu.

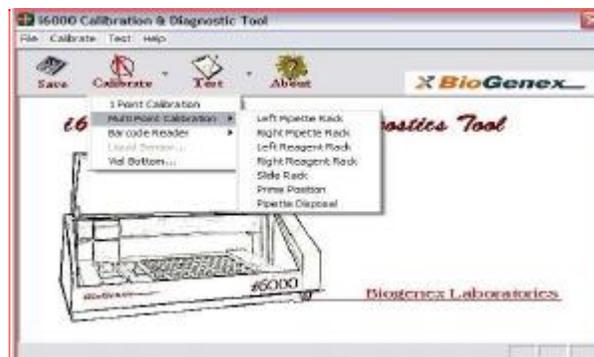


11.3 Left Pipette Rack

- From the 'Calibrate' menu, select the '**Multi Point Calibration**'

- Select **'Left Pipette Rack'** from the drop down list.
- The Z-head will move to position one of the left pipette tip box. The following screen will be displayed



- Make sure the Z-head is empty and place a pipette tip in position one in the left pipette box. Click 'OK' to continue.
- The Z-head will drop to pipette position one to pick up the pipette tip. The Calibrate Positions window is displayed on the screen

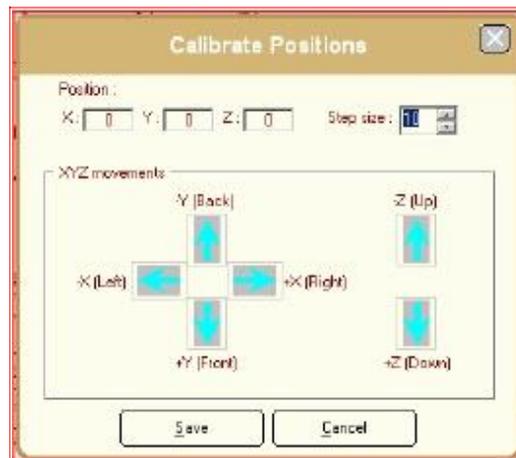

- Adjust the Z-head using the arrow keys as needed. Click **'Save'** to complete calibration.


To test the calibration for accuracy, use the diagnostics to pick up pipettes from positions 1, 8, 89, 96 and 45. To pass the test, pipette tips from all the positions should be picked up and be discarded in a single attempt.

11.4 Right Pipette Rack

- From the **'Calibrate'** menu, select the **'Multi Point Calibration'**
- Select **'Right Pipette Rack'** from the drop down list.
- The Z-head will move to position one of the right pipette tip box. The following screen will be displayed



- Make sure the Z-head is empty and place a pipette tip in position one in the right pipette tip box. Click 'OK' to continue.
- The Z-head will drop to pipette position one to pick up the pipette tip. The 'Calibrate Positions' window is displayed.

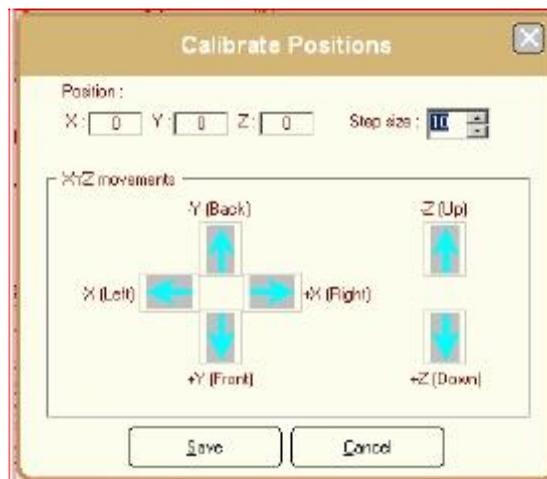


- Adjust the Z-head using the arrow keys as needed. Click '**Save**' to complete calibration.

To test the calibration for accuracy, use the diagnostics to pick up pipettes from positions 97, 104, 141, 185 and 192. To pass the test, pipette tips from all the positions should be picked up and be discarded in a single attempt.

11.5 Left Reagent Rack

- Place a vial in position one of the left reagent rack.
- Select 'Left Reagent Rack' from the 'Multi Point Calibration' drop down menu. The 'Calibrate Positions' window is displayed.

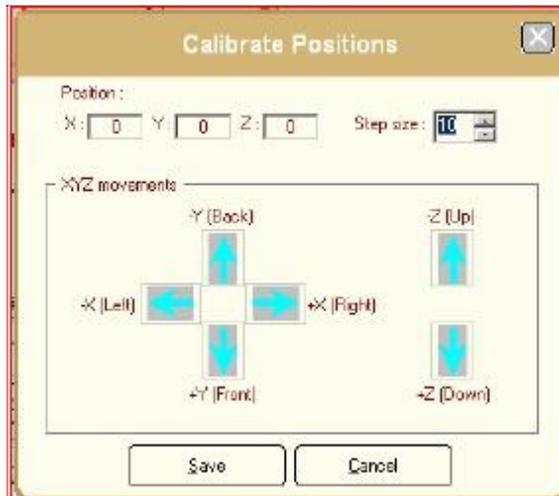

- Center the pipette tip adaptor over the reagent vial by using the X-Y axis control.
- Use the Z control to lower the pipette tip adaptor close to the vial to determine the center position. Click 'Save'

11.6 Right Reagent Rack


- Place a vial in position one of the right reagent rack.
- Select 'Right Reagent Rack' from the 'Multi Point Calibration' drop down menu. The 'Calibrate Positions' window is displayed.
- Center the pipette tip adaptor over the reagent vial by using the X-Y axis control.
- Use the Z control to lower the pipette tip adaptor close to the vial to determine the center position. Click 'Save'

11.7 Slide Rack

- Select 'Slide Rack' from the 'Multi Point Calibration' drop down menu.
- The Z-head will move towards slide position one. The following screen will be displayed.



- Make sure the Z-head is empty and click 'OK' to continue. The 'Calibrate Positions' window is displayed.

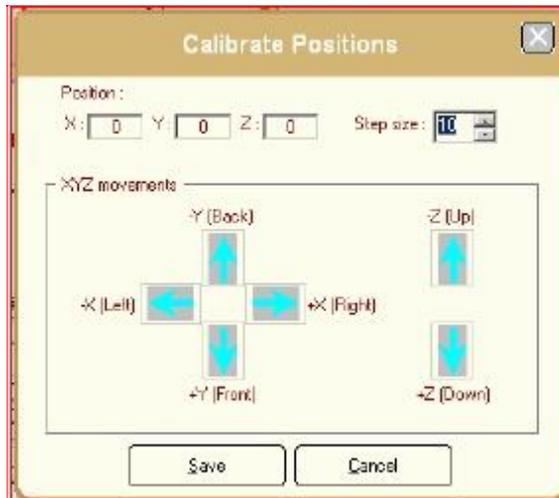

- Adjust the Z-head using the arrow key to ensure that the blow head is centered on the slide and located right after the label.
- Click '**Save**' to complete calibration.
To test the calibration for accuracy, use the diagnostics to blow and wash the slide. Make sure that the volume used covers the slide properly.

11.8 Prime Position

- Select Prime Position from the 'Multi Point Calibration' drop down menu.
- The Z-head will move forward to the front of the instrument.
- Adjust the Z-head position using the arrow key to ensure that the Z-head is far enough forward to allow the unit to prime the solutions without splashing the liquid on the slides. Save the settings.



11.9 Pipette Tip Disposal

- Select '**Pipette Disposal**' from the '**Multi Point Calibration**' drop down menu.
- The Z-head will move to a position above the pipette tip disposal area.
- Adjust the Z-head position using the arrow key so that the pipette tip can be dropped directly into the center of the disposal area. Save the settings.



11.10 Vial Bottom Calibration

Using a micro pipette, fill an OptiMiser™ Vial with 3000 µl of liquid and place the vial at position one. Be sure that the vial and rack are seated completely.

Select '**Test**' from the i6000™ Calibration & Diagnostic Window and then click on '**Diagnostics**' from the drop down list. The '**i6000 Diagnostics**' window is displayed.

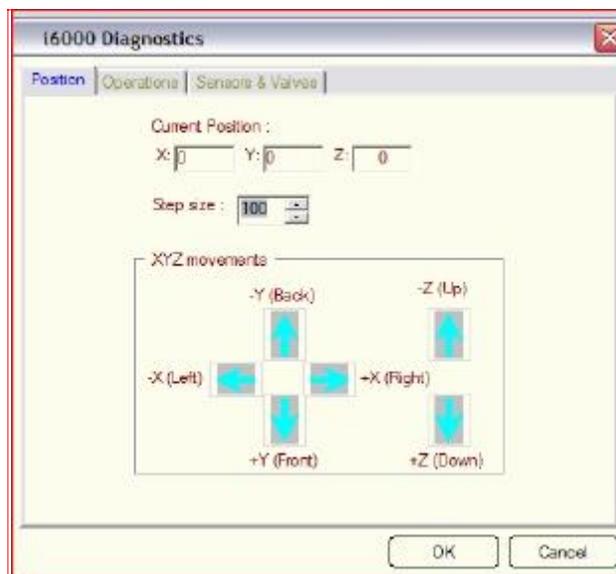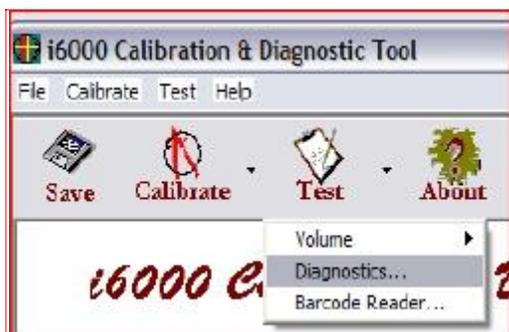


Click on the '**Operations**' and select '**Pickup Pipette**' from the options available.

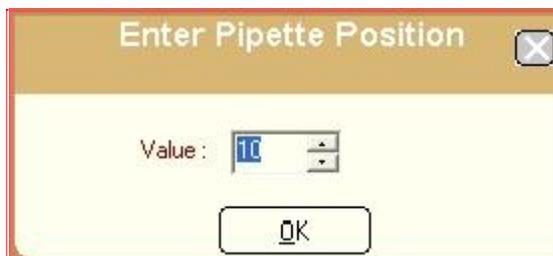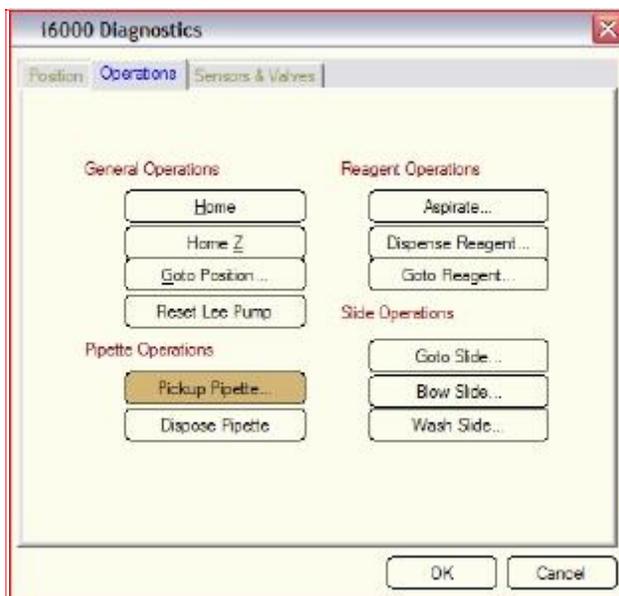

Enter a pipette position value and click '**OK**'. The Z-head will move forward to the pipette tray and pick up the pipette from the desired position. When the operation is done, close the '**Diagnostics**' window.

Open the 'i6000 Calibration and Diagnostic tool' window, and select '**Vial Bottom**' from the '**Calibrate**' menu.

- From the '**Calibrate**' menu, click on Vial Bottom from the drop down list.
- The Calibrate Reagent Vial Bottom dialog box will be displayed.

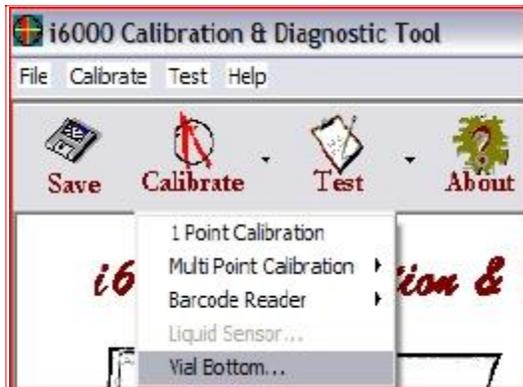

- The instrument will sense the liquid three times and display the value reading each time. After three attempts, the best sensing level is displayed





- Continue to click on **'Calibrate'** until a consistent reading is obtained three times in a row.
- Once a consistent reading is obtained, click on the done button to save the value in the registry.
- If a one-point calibration is performed, then the bottom vial calibration will need to be repeated.

11.11 Diagnostics

The i6000™ Infinity diagnostics feature is a powerful tool designed to assist you in checking the components for proper function. This program provides a method for you to test the motors, valves and pumps used in the system. The ability to test these devices can greatly reduce the time taken to identify failures.

- Open the 'i6000™ Calibration and Diagnostic tool' window, and select **'Diagnostics'** from the 'Test' menu.
- The Diagnostics screen will be displayed, providing you with complete control of the system from one convenient screen.

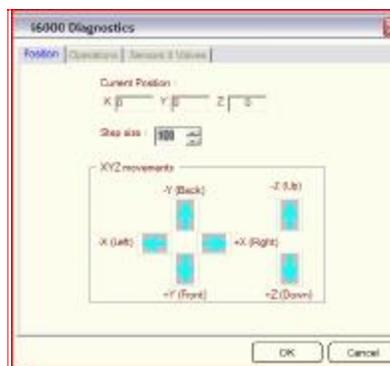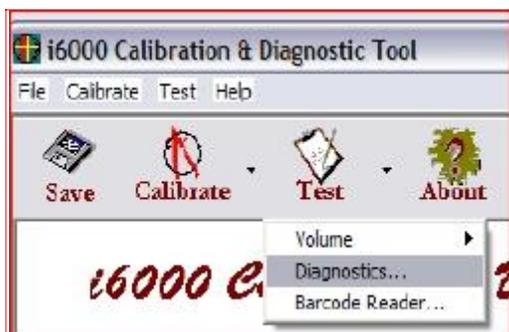


- The Diagnostic screen contains 3 main sections: Position, Operations and Sensors & Valves.

11.12 Position

- This section enables you to move the Z-head in any direction in the specified step size (1-1000).
- The step size determines how far the Z-head will move in the selected direction. Type in the step size then click on the arrow corresponding to the direction that you want to move.
- You can also use the 'Arrow Keys' on keyboard to move the Z-head.
- This section also displays the Current Position (expressed in number of steps from home) of the Z-head.

11.13 Operations

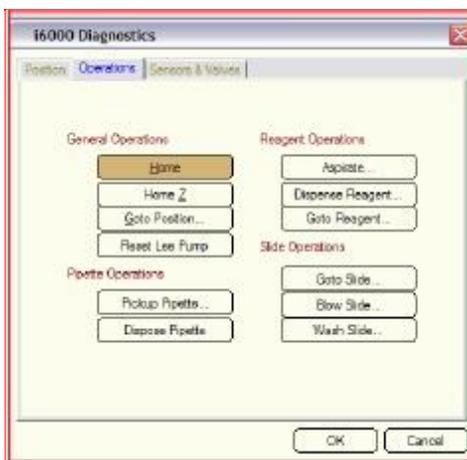

The following pre-programmed functions are available:

General Operations

- Home: Move Z-head to Home in all three directions X, Y, and Z.
- Home Z: Move the Z-head to home in the Z direction only.
- Go to Position: Allows you to enter the X, Y, and Z coordinates that you would like to move to.

- Reset Lee Pump: This is reserved for service personnel. It allows correct aspiration and dispensing of reagent.

Reagent Operations

- Aspirate: Move down in the current X and Y position until fluid is detected. Draw specified volume into pipette (1-1000).
- Dispense Reagent: In the current X, Y and Z position dispense specified volume (1-1000).
- Go to Reagent: Move to a specified reagent vial (1-60).


Pipette Operations

- Pickup Pipette: Pickup specified pipette tip. (1-192).
- Dispose Pipette: Move to the Disposal Position and eject tip from adapter.

Slide Operations

- Go to Slide: Move to a specified slide (1-60).
- Blow Slide: blow the slide
- Wash Slide: wash the slide

11.14 Sensors and Valves

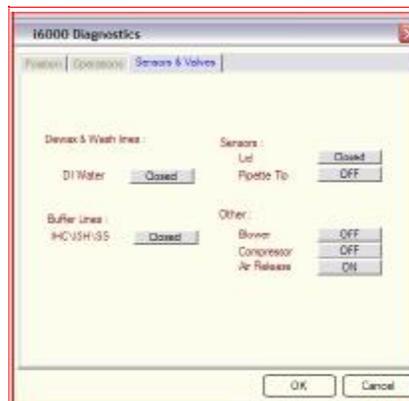

This section enables you to control the state (on or off) of the fluid delivery valves, blower, compressor and the air release valve. It also displays the status of the lid and the pipette tip sensor.

To test the fluid delivery valves you must turn the air release valve off and the compressor on and allow the system to pressurize.

Make sure the Z-head is moved over to the Prime Position. Once pressurized open the fluid delivery valve to be tested. Fluid should now be dispensed from the wash head.

12.0 Error Conditions

This section describes some of the common error messages that might appear while using the i6000™ Automated Staining System.


12.1 System Door Error

The instrument will not run with the door open. If the system door is open when starting the run, an error message will be displayed.



Close the system door and click **OK** to continue.

If at any time during the run, the system door is opened, the instrument will pause and an error message will be displayed.

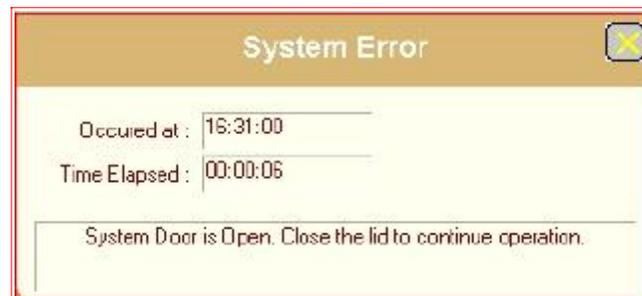

When the system door is closed the run will resume.

12.2 Pipette Errors

Pipette tip errors are divided into four basic groups: Tip Errors at Box, Tip Errors at Reagent, Tip Errors at Slide, and Tip Errors Post-Dispense. Each group provides a specific error message and requires a slightly different response.

12.2.1 Pipette Tip Errors at Box

- If the Z-head tip adapter fails to pick up a pipette tip from the box holder, or the pipette tip is missing, the system will re-attempt to pick up a tip.
- The Z-head goes to its home position, returns to the pipette tip it just failed to pick up, and tries to pick it up again.
- If it successfully picks up the tip on this second try at the same position, the staining run continues. If it fails to pick up the tip at the same position, it will attempt to pick up a tip in the next three positions.

- After four consecutive failures, an audible alarm will sound and the following message will be displayed.



- Replace the box of pipette tips in the unit with a new full box.
- Make sure the tip box is well seated and that there is no obvious obstruction to the Z-head movement. Select Resume continuing.
- When Resume is selected, the Z-head will attempt to pick up the first tip in the new box. If it is successful, the staining run continues. If not, the cycle begins again.
- The Z-head will make four attempts to pick up a tip from the new box of pipette tips. If unsuccessful, the above error message will be displayed again.
- Click '**Cancel the Run**' to abort the run.

12.2.2 Pipette Tip Errors at Reagent

If the Z-head tip adapter drops a pipette tip after drawing reagent, but before moving to the slides, the Z-head will stop, an alarm sets off with the following message displayed on the screen.



A 10 sec. timer will start ticking. Remove the pipette tip and click '**OK**' to continue the run. If it takes more than 10 sec. to remove the tip, a 'Critical Error' message will be displayed on the screen.



Make sure the reagent vial and the vial rack are positioned properly. Follow the instructions provided in the error message, and click **'Resume'**.

The Z-head will pick up a new pipette tip and move to the reagent vial to draw the reagent again and continue the run.

12.2.3 Pipette Tip Errors at Slide

If the Z-head tip adapter drops a pipette tip after moving to the slides but before completely dispensing the drawn reagent, the Z-head will stop, and an alarm sets off with the following message displayed on the screen.



If the error is not cleared within 10 sec. a 'Critical Error' message will be displayed on the screen.




Make sure the slide racks are positioned properly and there is no other obstruction to the free movement of the Z-head. Follow the instructions provided in the error message and click 'Resume'.

The Z-head will pick up a new pipette tip, to dispense the reagent on each of the affected slides listed in the error message.


12.2.4 Pipette Tip Errors Post Dispense

If the Z-head tip adapter drops a pipette tip after completely dispensing the reagent, but before the tip is ejected, the Z-head will stop.

An alarm sets off with the following message displayed on the screen.

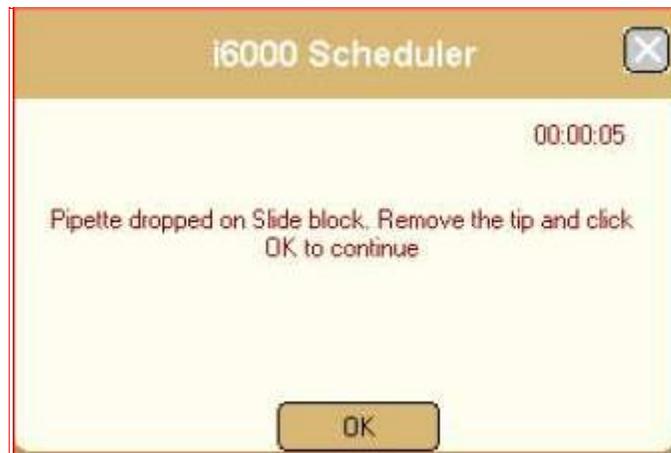

If the error is not cleared within 10 sec. a 'Critical Error' message will be displayed.

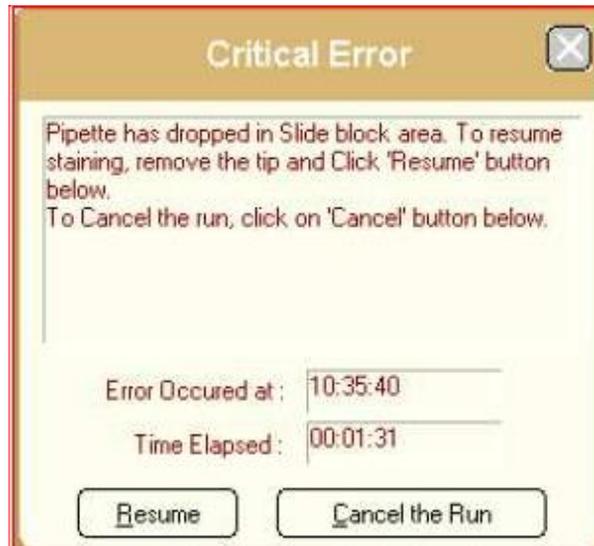


Remove the dropped tip and select Resume to continue the run.

12.3 Reagent Errors

Reagent Vial Empty or Missing Errors

If the reagent vial is either empty or does not contain the required amount of reagent, the liquid level sensor will detect and sets off an alarm with the following message displayed on the screen.

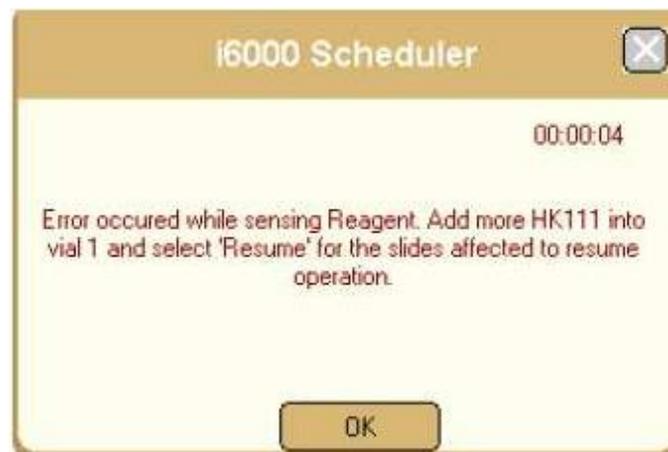

All the slides that are affected are highlighted in a different color (selected in the color scheme), and the staining is stopped at this point for the affected slides.

These slides will then be buffered periodically depending upon the time set in global settings for Buffer interval during run. To resume the staining, right click on the slide and select '**Resume staining**'.

12.3.1 Insufficient Reagent Vials Error

If a run requires more than 60 reagents, including primary antibodies, links, labels, substrates, and counter-stains, clicking on **'Start'** will display the following message.



- Click **'OK'** to cancel the run.
- You can start over again or in the case of an open run, open the previous stored staining run and edit it.
- Edit/ Delete slides and/or protocols to bring down the total number of reagents required to 60 or less.

12.4 System Pressurization Errors

If the system is not pressurized within five minutes, the **'i6000 Scheduler'** dialog box will be displayed.



Ensure that all carboy caps are securely tightened and click on the **'Retry'** button.

If pressurization continues to fail, click on the **'Cancel the Run'** button to exit. Check all carboy connections and secure them before restarting the run.

12.5 Substrate and Label Error

When the substrate selected does not match the label selected, an error message will appear when you try to save the protocol.

Press the **Enter** key, or click '**OK**' to return to the **Edit** screen and make the correction.



Appendix-1

System Specifications

- Dimensions:
40.5" W x 24" D x 18.5" H (1028 mm W x 610 mm D x 470 mm H)
- Weight:
130lbs
- Electrical Specifications:
115/230VAC 50/60Hz
- Power Consumption:
300 Watts
- Normal Operating Temperature:
18-28.5 °C (64 – 80.5 °F)
- Slide Capacity:
Total Capacity: 60 Slides, 5 racks (12 slides per rack)
- Reagent Capacity:
Total Capacity: 60 Reagent vials, 2 racks (30 reagents vials per rack)
- Reagent Dispensing Volume:
100 µl, 150 µl, 200 µl, 250 µl, 300 µl, 350 µl, 400 µl, 500 µl, 600 µl, 700 µl, 800 µl and 900 µl
- Reagent Dispensing Accuracy: ± 5 %
- Operating system: Windows® XP
- Buffer Containers: 10 Liters Each
- Waste Containers: 20 Liters Each

Performance Characteristics

The pattern and intensity of staining is highly consistent and reproducible when slides from the same run using identical specimen and protocols are compared. Similarly, slides using the same specimen and protocols stained on different instruments yielded highly consistent and reproducible results. In comparison to manual results, the instrument staining results are identical, if not indistinguishable. Staining specifications were consistently achieved in STAT, continuous or batch run formats.

Appendix-2**Service Plan**

For a period of one year from the date of installation ("Service Period"), BioGenex Laboratories, Inc. ("BioGenex") will provide the following services for maintenance and repair of the i6000™ Automated Staining System ("System") purchased, leased or rented by the Customer ("Services"):

1. Preventative Maintenance Service
2. Emergency Repair Service
3. Software Upgrades

Customer Responsibilities

Customer must (1) maintain the installation site and its environment in a condition suitable for operation of the System, and (2) maintain the System under the daily/monthly routine maintenance schedule as provided in the System Operating Manual.

Service Charges

BioGenex may charge Customers its standard service rates or may decline to provide the services required to correct a malfunction caused by Customer's failure to fulfill its responsibilities hereunder, the failure of anyone other than BioGenex or its service contractor to comply with its written instructions or recommendations, the combination of the System with an incompatible third party product, the alteration or improper storage, handling, use or maintenance of any part of the System by anyone other than BioGenex or its service contractor, any factor external to the System or beyond BioGenex reasonable control.

Additional Terms

Services will be provided during regular business hours. BioGenex may subcontract with service contractors any of the service obligations to Customer. No such subcontract will release BioGenex from its obligations to Customer. In lieu of repair, BioGenex may elect to replace the System or any part thereof, or refund the portion of the Purchase Price paid by Customer. Replacement parts will be provided on an exchange (refurbished) or new part basis, at BioGenex's option.

Labor to install replacement parts is included. Replaced parts become BioGenex's property. Components of the System covered by a third party warranty are subject to the limitations contained in that warranty. The Services are due only if the System remains at the Installation Site. Rigging or facility services, accessory, consumable and supply items are charged separately. The Services provided under this Service Plan are performed under the terms and conditions of the Standard System Purchase and License Agreement or Lease and Reagent Purchase Agreement between BioGenex and Customer, as applicable, which are incorporated herein by reference.

Warranty Disclaimer

The System, including all hardware, software, and parts are provided AS IS. BIOGENEX SPECIFICALLY DISCLAIMS ALL EXPRESS, STATUTORY OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT OF THIRD PARTY RIGHTS. Where disclaimer of implied warranties is prohibited by law, such warranties expire 90 days after the date of delivery of the System. Customer acknowledges and accepts that operation of the Software may not be uninterrupted or error free and that any warranty for Year 2000 compliance is made only to the extent provided in the software licenses contained in this Operating Manual.

Appendix-3

Software License (BioGenex Software)

IMPORTANT- READ CAREFULLY BEFORE USING THE i6000™ AUTOMATED STAINING SYSTEM WHICH CONTAINS MICROSOFT® SOFTWARE.

This is a license and not an agreement for sale. You obtain no rights other than those granted you under this BioGenex Laboratories Inc., ("BioGenex") Software License ("License"). By using the i6000™ Automated Staining System on which BioGenex software program(s) have been preinstalled, you are agreeing to be bound by the terms of this license.

1. LICENSE.

The term "Software" means the programs that are installed in and supplied with a BioGenex computerized system ("System"), such as an i6000™ Automated Cell System, and subject to this License, but excluding the Microsoft software, which is subject to a separate license. BioGenex grants you a non-exclusive license to use the Software only in the System in which such Software is initially installed, and only for so long as you are using such System at the Site approved by BioGenex. BioGenex will provide one copy of the Software with the initial installation of the System, some of which may be preinstalled on the System's hard drive. You may make one backup copy, which is subject to the terms of this License and which must include the copyright notice(s) and any legend(s). All copies of the Software remain the property of BioGenex or its affiliates or suppliers. No title, ownership, or other interest in the Software is transferred to you, except as expressly provided herein. You will not modify, reverse engineer, decompile, disassemble, or derive source code from, the Software. In countries which permit reverse engineering notwithstanding the previous sentence, you agree not to conduct such reverse engineering until you have notified BioGenex in writing and given BioGenex a reasonable time to provide such information as you are entitled by law.

Unless BioGenex specifies otherwise, any modifications to, or new releases or versions of, the Software will be provided subject to this License at applicable charges.

BioGenex will grant a license to the party to whom you transfer possession of the Software ("Subsequent Licensee") only when all of the following terms are met:

(1) BioGenex approves such transfer in writing in its sole discretion; **(2)** the Software, including all copies, and the System are concurrently transferred to the Subsequent Licensee; and **(3)** the Subsequent Licensee agrees in writing to be bound by the terms of this Software License with

respect to the Software transferred. Transfer of the Software terminates your license from BioGenex.

2. INTELLECTUAL PROPERTY.

The Software is the confidential trade secret of BioGenex. You agree to take precautions to prevent any disclosure or use of Software (except as specifically permitted herein) consistent with precautions used to protect your own confidential trade secrets, but in no event less than reasonable care.

3. WARRANTY.

If the Software is warranted, the terms and conditions of such warranty is specified in the agreement under which such Software was originally acquired from BioGenex or a BioGenex Authorized Distributor ("Distributor"). Otherwise, the Software is provided "AS IS." The operation of the Software may not be uninterrupted or error free. SUCH WARRANTIES REPLACE ALL OTHER WARRANTIES, EXPRESS, STATUTORY OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT OF THIRD PARTY RIGHTS. Support services for Software which is not under warranty may be available from BioGenex or a Distributor for a charge.

4. LIMITATION OF LIABILITY.

THE MAXIMUM AGGREGATE LIABILITY OF BIOGENEX AND ITS SUPPLIERS, DEVELOPERS, SUBCONTRACTORS, AND AGENTS FOR ANY CLAIMS ARISING OUT OF OR RELATED TO THIS AGREEMENT, WHETHER IN CONTRACT, TORT, PRODUCT LIABILITY OR OTHERWISE, SHALL NOT EXCEED THREE THOUSAND DOLLARS.

Under no circumstances will BioGenex be liable for losses or damages resulting from (1) third party claims against you or claims by you based on third party claims (other than those referred to in item 1 in the first paragraph of this Section), (2) loss of stored, transmitted or recorded data, (3) consequential damages (including lost profits or savings), or incidental, indirect or punitive damages, even if BioGenex or a Distributor is informed of their possibility, or (4) inability to perform under this License, in whole or in part, caused by acts of God or any other event beyond BioGenex's control. This Section sets forth the maximum collective responsibility of BioGenex, its suppliers, developers, subcontractors and agents, and all such parties are intended beneficiaries of this Section.

5. U.S. GOVERNMENT RESTRICTED RIGHTS.

If you are a U.S. Government user, the Software and documentation are commercial computer software developed exclusively at private expense and are subject to the terms of this License pursuant to 48 C.F.R. 12.212 or 48 C.F.R. 227.7202. The manufacturer is BioGenex Laboratories Inc., 48810 Kato Road, Suite 100E & 200E, Fremont, CA 94538 Tel : +1 (800) 421-4149 Fax : +1 (510) 824-1490 E-mail: support@biogenex.com.

6. GENERAL.

You may terminate your license at any time by destroying all your copies of the Software or as otherwise described in this License. BioGenex may terminate your license if you do not comply with the terms of this License. Upon such termination, you will destroy all your copies of the Software. Any attempt to sublicense, rent, lease, assign or, except as expressly provided herein, to transfer any copy of the Software is void. This License is governed by the laws of the State of California without regard to its conflicts of law rules.

Appendix-4

Software License Agreement (Microsoft Software)

IMPORTANT-READ CAREFULLY: This End-User License Agreement ("EULA") is a legal agreement between you (either an individual or a single entity) and the manufacturer ("Manufacturer") of the special purpose computing device ("SYSTEM") you acquired which includes certain Microsoft software product(s) installed on the SYSTEM ("SOFTWARE PRODUCT" or "SOFTWARE"). The SOFTWARE includes computer software, the associated media, any printed materials, and any "online" or electronic documentation. By installing, copying or otherwise using the SOFTWARE, you agree to be bound by the terms of this EULA. If you do not agree to the terms of this EULA, Manufacturer and Microsoft Corporation ("Microsoft") are unwilling to license the SOFTWARE to you. In such event, you may not use or copy the SOFTWARE, and you should promptly contact Manufacturer for instructions on return of the unused product(s) for a refund.

SOFTWARE LICENSE

The SOFTWARE is protected by copyright laws and international copyright treaties, as well as other intellectual property laws and treaties. The SOFTWARE is licensed, not sold.

1. GRANT OF LICENSE. This EULA grants you the following rights:

Software: You may use the SOFTWARE as installed on the SYSTEM.

Application Sharing: The SOFTWARE may contain technology that enables applications to be shared between two or more SYSTEMS, even if an application is installed on only one of the SYSTEMS. You should consult your application license agreement or contact the application licensor to determine whether sharing the application is permitted by its licensor.

Storage/Network Use: If the SOFTWARE PRODUCT is installed on the SYSTEM over an internal network from a server, you must acquire and dedicate a license for the SOFTWARE PRODUCT for each SYSTEM on which the SOFTWARE PRODUCT is used or to which it is distributed. A license for the SOFTWARE PRODUCT may not be shared or used concurrently on different SYSTEMS.

Back-up Copy: If Manufacturer has not included a copy of the SOFTWARE on separate media with the SYSTEM, you may make a single copy of the SOFTWARE for use solely for archival purposes with the SYSTEM.

2. DESCRIPTION OF OTHER RIGHTS & LIMITATIONS.

Limitations on Reverse Engineering, Decompilation and Disassembly: You may not reverse engineer, decompile, or disassemble the SOFTWARE, except and only to the extent that such activity is expressly permitted by applicable law notwithstanding this limitation.

Single SYSTEM: The SOFTWARE is licensed with the SYSTEM as a single integrated product. The SOFTWARE may only be used with the SYSTEM.

Rental: You may not rent or lease the SOFTWARE.

Software Transfer: You may permanently transfer all of your rights under this EULA only as part of a sale or transfer of the SYSTEM, provided you retain no copies, you transfer all of the SOFTWARE (including all component parts, the media, any upgrades or backup copies, and this EULA, and if applicable, the Certificate(s) of Authenticity), **and** the recipient agrees to the terms of

this EULA. If the SOFTWARE is an upgrade, any transfer must include all prior versions of the SOFTWARE.

Termination: Without prejudice to any other rights, Manufacturer or Microsoft may terminate this EULA if you fail to comply with the terms and Conditions of this EULA. In such event, you must destroy all copies of the SOFTWARE and all of its component parts.

Single EULA: The package for the SOFTWARE may contain multiple versions of this EULA, such as multiple translations and/or multiple media versions (e.g. in the user documentation and in the software). In this case, you are only licensed to use one (1) copy of the SOFTWARE PRODUCT.

3. UPGRADES AND RECOVERY MEDIA

If the SOFTWARE is provided by Manufacturer on media separate from the SYSTEM and is labeled "For Upgrade Purposes Only" ("Upgrade SOFTWARE"), you may install one copy of the Upgrade SOFTWARE onto the SYSTEM as a replacement copy for the SOFTWARE originally installed on the SYSTEM and use it in accordance with Section 1 of this EULA. You may also install additional copies of the Upgrade SOFTWARE as replacement copies onto additional SYSTEMS which are the same brand and model as the SYSTEM and contain a duly licensed copy of the same version and language release of the SOFTWARE ("ADDITIONAL SYSTEMS"), provided that (1) Manufacturer has supplied a corresponding serialized sticker for each additional copy of the Upgrade SOFTWARE, and (2) you affix a serialized sticker per Manufacturer's instructions for each unit of Upgrade SOFTWARE you install. If the SOFTWARE is provided by Manufacturer on separate media and labeled as "Recovery Media," you may not make a copy of the SOFTWARE as described in Section 1 for archival purposes. Instead, you may use the Recovery Media solely to restore or reinstall the same version and language release of the SOFTWARE as originally installed on the SYSTEM and thereafter use the SOFTWARE as restored or reinstalled in accordance with Section 1 of this EULA. A single unit of Recovery Media may be used by you to restore or reinstall the SOFTWARE on ADDITIONAL SYSTEMS.

4. COPYRIGHT

All title and copyrights in and to the SOFTWARE (including but not limited to any images, photographs, animations, video, audio, music, text and "applets," incorporated into the SOFTWARE), the accompanying printed materials, and any copies of the SOFTWARE, are owned by Microsoft or its suppliers. You may not copy the printed materials accompanying the SOFTWARE. All rights not specifically granted under this EULA are reserved by Microsoft.

5. PRODUCT SUPPORT

Product support for the SOFTWARE is not provided by Microsoft or its subsidiaries. For product support please refer to Manufacturer's support number provided in the documentation for the SYSTEM. Should you have any questions concerning this EULA, or if you desire to contact Manufacturer for any other reason, please refer to the address provided in the documentation for the SYSTEM.

6. LIMITED WARRANTY

Limited Warranty: Manufacturer warrants that the SOFTWARE will perform substantially in accordance with the accompanying written materials for a period of ninety (90) days from the date of receipt. Any implied warranties on the SOFTWARE are limited to ninety (90) days. Some states/jurisdictions do not allow limitations on duration of an implied warranty, so the above limitation may not apply to you.

Customer Remedies: Manufacturer's and its suppliers' entire liability and your exclusive remedy shall be, at Manufacturer's option, either (a) return of the price paid, or (b) repair or replacement of the SOFTWARE that does not meet the above Limited Warranty and which is returned to Manufacturer with a copy of your receipt. This Limited Warranty is void if failure of the SOFTWARE has resulted from accident, abuse, or misapplication. Any replacement SOFTWARE will be warranted for the remainder of the original warranty period or thirty (30) days, whichever is longer.

No Other Warranties. EXCEPT AS EXPRESSLY PROVIDED IN THE LIMITED WARRANTY SECTION ABOVE, THE SOFTWARE IS PROVIDED TO THE END USER "AS IS" WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, WARRANTIES OF NON-INFRINGEMENT, MERCHANTABILITY, AND/OR FITNESS FOR A PARTICULAR PURPOSE. THE ENTIRE RISK OF THE QUALITY AND PERFORMANCE OF THE SOFTWARE IS WITH YOU.

No Liability for Consequential Damages. MANUFACTURER'S SUPPLIERS SHALL NOT BE HELD TO ANY LIABILITY FOR ANY DAMAGES SUFFERED OR INCURRED BY THE END USER (INCLUDING, BUT NOT LIMITED TO, GENERAL, SPECIAL, CONSEQUENTIAL OR INCIDENTAL DAMAGES INCLUDING DAMAGES FOR LOSS OF BUSINESS PROFITS, BUSINESS INTERRUPTION, LOSS OF BUSINESS INFORMATION AND THE LIKE), ARISING FROM OR IN CONNECTION WITH THE DELIVERY, USE OR PERFORMANCE OF THE SOFTWARE.

7. NOTE ON JAVA SUPPORT.

THE SOFTWARE PRODUCT MAY CONTAIN SUPPORT FOR PROGRAMS WRITTEN IN JAVA. JAVA TECHNOLOGY IS NOT FAULT TOLERANT AND IS NOT DESIGNED, MANUFACTURED, OR INTENDED FOR USE OR RESALE AS ON-LINE CONTROL EQUIPMENT IN HAZARDOUS ENVIRONMENTS REQUIRING FAIL-SAFE PERFORMANCE, SUCH AS IN THE OPERATION OF NUCLEAR FACILITIES, AIRCRAFT NAVIGATION OR COMMUNICATION SYSTEMS, AIR TRAFFIC CONTROL, DIRECT LIFE SUPPORT MACHINES, OR WEAPONS SYSTEMS, IN WHICH THE FAILURE OF JAVA TECHNOLOGY COULD LEAD DIRECTLY TO DEATH, PERSONAL INJURY, OR SEVERE PHYSICAL OR ENVIRONMENTAL DAMAGE.

U.S. GOVERNMENT RESTRICTED RIGHTS

The SOFTWARE and documentation are provided with RESTRICTED RIGHTS. Use, duplication or disclosure by the Government is subject to restrictions as set forth in subparagraph (c)(1)(ii) of the Rights in Technical Data and Computer Software clause at DFARS 252.227-7013 or subparagraphs (c)(1) and (2) of the Commercial Computer

Software - Restricted Rights at 48 CFR 52.227-19, as applicable. Manufacturer is Microsoft Corporation/One Microsoft Way/Redmond, WA 98052-6399.

If you acquired this EULA in the United States, this EULA is governed by the laws of the State of Washington. If you acquired this EULA in Canada, this EULA is governed by the laws of the Province of Ontario, Canada. Each of the parties hereto irrevocably attorns to the jurisdiction of the courts of the Province of Ontario and further agrees to commence any litigation, which may arise hereunder in the courts located in the Judicial District of York, Province of Ontario. If this EULA was acquired outside the United States, then local law may apply. Should you have any questions concerning this EULA, please contact the Manufacturer.

Appendix-5

Daily /Monthly Cleaning

Daily Cleaning

The i6000 Automated Staining System requires minor cleaning and sanitation to ensure consistent staining, eliminate potential contamination, and guarantee the longevity of the instrument. The following procedure should be performed after the last staining run of the day:

1. Remove the Pipette Tip Waste Tray. Discard tips and rinse tray with tap water.
2. Remove the tip boxes, reagent rack and slide racks from the instrument. Rinse the reagent rack and slide racks with tap water.
3. Disconnect the buffer carboy from the connector panel of the instrument.
4. Connect the DI Water Carboy to buffer quick coupling (White color) on the liquid panel of the instrument.
5. Prime the system with buffer. The i6000 Automated Staining System will rinse the buffer line with the DI water. Repeat three times.
6. Now disconnect the DI Water carboy from the buffer quick coupling and connect it to the DI Water quick coupling (Green color) on the liquid panel.
7. Prime the system with DI Water. The i6000 Automated Staining System will rinse the DI Water line with the DI water. Repeat three times.
8. After the priming cycle is finished, remove the rinse tubing and reconnect each carboy to its appropriate inlet port.
9. Wipe the working space of the instrument with a soft damp cloth. Dry the interior by wiping all surfaces with a clean soft cloth.

Monthly Cleaning

To prevent the growth of bacteria or mold in the tubing, carboys, and other areas, the following procedure should be performed monthly. Wear gloves and eye protection while performing this procedure:

1. Remove the tip boxes, reagent rack and slide racks from the machine. Rinse the reagent rack and slide racks with tap water.
2. Rinse the tip adapter, tip ejector and the blow head with steady stream of DI water from a rinse bottle. Remove any salt accumulation.
3. Remove the carboys and caps from the system including the Waste carboy. Empty both the Buffer and DI Water carboys and rinse with DI Water. Fill each carboy with 1 Gallon of a 1% bleach solution. (30 - 40ml of bleach to 4 liters of DI Water). Secure the carboys with caps and roll the bleach solution thoroughly around the carboys.
4. Empty the bleach solution from each of the carboys and rinse with DI Water. Empty the carboys.
5. Store the Buffer and DI Water carboys upside down to let them dry.
6. In the DI Water carboy, prepare fresh cleaning solution by diluting 150 ml of Professional Lysol Deodorizing cleaner (active ingredient: alkyl dimethyl benzyl ammonium chlorides) to 3600 ml of distilled water. In Countries, where Professional Lysol is not available or is not

allowed to be sold, use 4% Sodium Hypochlorite. Dilute 50ml of Sodium Hypochlorite in 950ml of DI Water.

7. Connect the Waste Carboy to the instrument.
8. Connect the DI Water Carboy filled with diluted solution to buffer quick coupling (White color) on the liquid panel of the instrument.
9. In prime solutions window, select the buffer line. The i6000 Automated Staining System will rinse the buffer line. Repeat three times.
10. Now disconnect the DI Water carboy filled with the diluted solution to the buffer quick coupling and connect it to the DI Water quick coupling (Green color) on the liquid panel.
11. In prime solutions window, select the DI Water line. The i6000 Automated Staining System will rinse the DI Water line. Repeat three times.
12. Wait for 10 minutes after the last prime is done.
13. Empty the DI Water carboy to discard any residual solution and rinse the carboy with DI Water.
14. Fill the DI Water carboy with DI and rinse both the lines by priming 3 times.
15. Empty the DI Water carboy and store it upside down and let it dry.
16. Clean the Drain Tray with powdered cleanser (Comet with bleach) and a soft plastic brush. Rinse and dry. Check that the waste drains properly. Adjust the drain tubing as necessary.
17. Clean the working space table with a soft damp cloth.
18. Dry the interior of the instrument by wiping all surfaces with a clean soft cloth.
19. Leave the lid open and let the instrument air dry overnight.

Appendix-6

Instrument Packaging Procedure for i6000™ Diagnostics/ Infinity Automated Staining System

1.0 Purpose:

1.1 This procedure describes the packing procedure for i6000™ Diagnostics/ Infinity shipment.

2.0 Scope:

2.1 This procedure applicable to the Systems Manufacturing Department.

3.0 Definitions:

3.1 QC: Quality Control

4.0 Responsibility:

4.1 It is the responsibility of Systems Manufacturing to oversee and enforce this procedure.

5.0 Related Documents:

5.1 Instrument Packaging Procedure Check list for Instrument Manufacturing Area doc no. 951-3320.1

5.2 Packing List for Accessories Box - i6000™ System (Doc. #951-3323.1)


6.0 Packaging procedure for shipping an instrument without accessory kits follows:

6.1 Two people are required when handling the instruments.

6.2 Inspect instrument to be sure it is undamaged and clean. All reagents, pipette tips, tip tray and slide carriers should be removed.

6.3 Place the instrument into a vapor bag (shrink-wrap is acceptable) and seal on the top.

6.4 Place the Cap (PN: 6100-20005) on the flat floor and place Fab, Bottom Buffer Left (PN: 6100-20004) and Fab, Bottom Buffer Right (PN: 6100-20003) inside the Cap carton (PN: 6100-20005) as shown in the photo1.

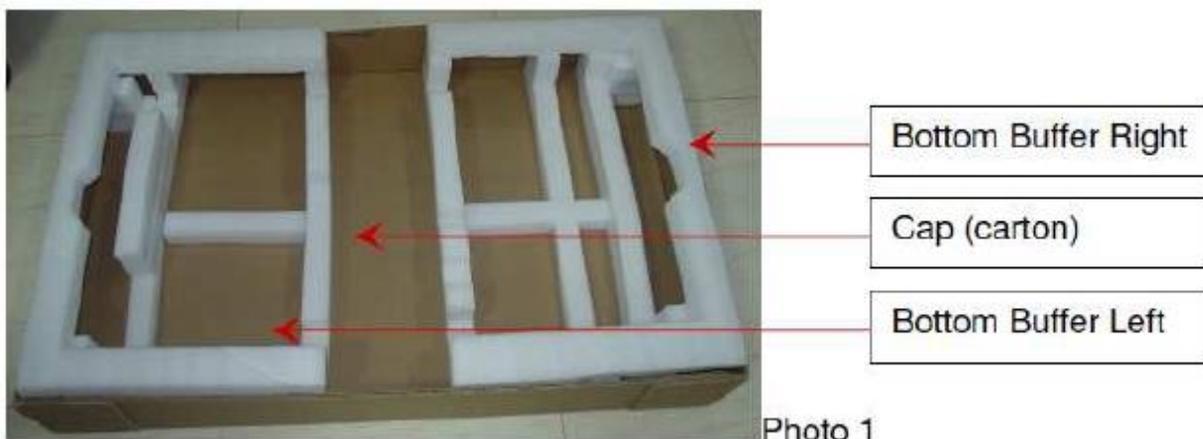

6.5 Place the i6000™ Autostainer instrument on the both the buffers and ensure that instrument has properly seated on the buffers. Place the Sleeve carton (PN: 6100-20005) from the top of the instrument as shown in the photo 2.



Instrument

Sleeve

Photo 2

6.6 Place the Fab, Buffer Right (PN: 6100-20001) and Fab, Buffer Left (PN: 6100-2002) on the instrument as shown per the photo 3.



Top Buffer Right

Top Buffer Left

Photo 3

6.7 Place the Top Cap on the sleeve as shown in the photo 4.




Photo 4

6.8 Wrap it with the nylon band strip using edge protectors. Push this setup into the Wooden Crate as shown in photo 5



Photo 5

6.9 Close the Crate with the front door and lock latches. Strap the door with metal banding material and clips in no less than 2 locations. The bands should run from end to end parallel to the long side of the pallet (see reference drawing). This is important to make sure the bands are not broken by pallet jacks or forklifts during shipping. Bands must be tight enough to prevent the instrument box from moving on the pallet.

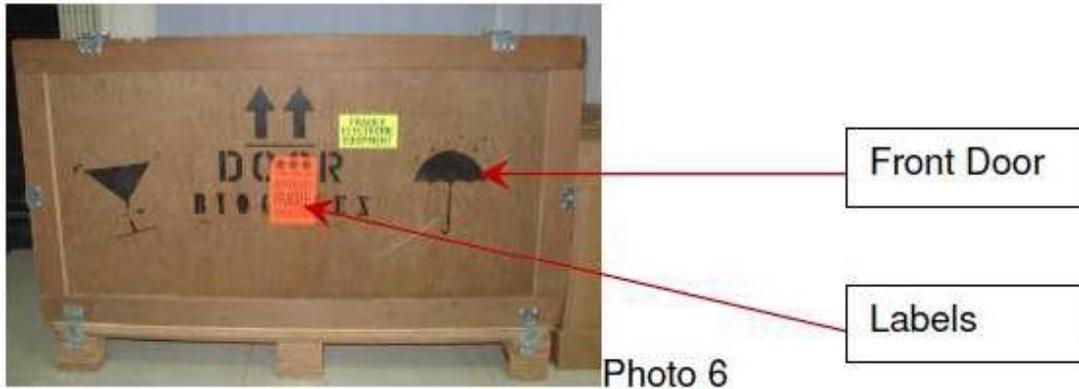

Apply appropriate container marking and labels as required by the Package Marking Procedure based on the shipments destination and mode of transportation.

7.0 The following instructions apply if the instrument is to be shipped with an accessory kit:

7.1 Two people are required when handling a loaded accessory kit box.

7.2 Verify the accessories in the Box per Doc.#951-3323.1

7.3 Include a copy of the Packing list 951-3320.1 in the Accessory kit box.

7.4 Seal the accessory kit box on top and bottom in the same manner as the instrument box.

7.5 Accessory box should be placed on the Wooden crate of the instrument and strap it to the instrument box at 2 locations from front to rear using the same approved metal banding materials as for the instrument box.

7.6 Apply labels shipping labels on the crate as shown in the photo6.

Appendix-7

Pre-Installation Checklist for i6000™ Systems

Instructions to the Sales Team

This document helps the Customer in understanding the requirements for installing i6000™ instrument in their laboratory. Being prepared with an understanding and compliance to these requirements will help the installation team to perform efficient and hassle-free installation and training at the Customer's site.

Please help the Customer fill out this checklist and check the boxes to show that the topic was discussed with the Customer and approved by them. Please complete this document with clear and legible notes of any special instructions or confirmation to the specifications. Please fax or send the completed document to Customer Support for approval, and leave a completed copy with the Customer.

Scheduling the installation and training cannot take place if this document is not approved by all concerned.

BioGenex Representative: _____

Name of the Institution: _____

Institutional Contact:

Name: _____

Title: _____

Telephone: _____

Fax: _____

Email: _____

1. Delivery

- The Customer shall inspect all shipments upon delivery, for external damage. If there is evidence of external damage, BioGenex Customer Support shall be notified immediately.
- Dock working hours and access restrictions, if any.
- It may require at least two people to take the instrument out of the carton and carry it to its final install site.
- Appropriate cart for transport of accessories.
- Storage space for the empty carton/packaging material.

Recommendations / Customer comments:

2. Environmental Requirements

a. Space requirements

Minimum Table/Platform Size:

- Length = 72" (includes Instrument, Computer and Monitor)
- Width = 30"
- Sufficient to support 75kg
- Storage space to install the bulk reagent carboys beneath the table/platform
- Clearance of 6" is required from the rear wall.

Instrument Height = 20" (Lid Closed); 33" (Lid Open)

Recommendations / Customer comments:

b. Power requirement

- 230V, 50Hz, 15 Amps grounded wall socket for 1KVA Online UPS System that supports the Instrument, Computer and Monitor.
- UPS System shall be provided by the customer. The output of the UPS System shall have 3 sockets of 5A each.
- Additional wall socket for DeskJet Printer is required as the UPS rating do not include this printer.

Recommendations / Customer comments:

c. Room temperature requirement

- Ambient temperature range : 18°C - 28.5 °C

Recommendations / Customer comments:

d. DI Water

- The instrument requires DI water to be filled in the bulk reagent bottle.

Recommendations / Customer comments:

e. Waste Disposal

- The instrument has a built-in waste storage container and does not require connecting to a running drain. Local regulations with respect to waste storage and disposal should be observed by the customer.

Recommendations / Customer comments:

3. Training information

- Key Operator/s to be trained: _____

- Is/Are the Key Operator/s competent with using MS WINDOWS Operating system?
- Please have the expected Key Operator/s answer the following questions
 1. Describe the use of the right mouse button
 2. Use the "Start" button to start Windows Explorer
 3. Arrange the desktop alphabetically
 4. Use the Clock to set date and time
 5. Open "Windows Explorer" and open a directory on the "C:" drive.
 6. Copy a file from one directory to another, using "drag and drop"
 7. Use "My Computer" to select a printer and display its properties.
 8. Open the "Recycle Bin".

- If the Key Operator/s cannot answer six of the eight questions, formal training on MS WINDOWS must be accomplished by the customer prior to installation and training with the i6000™ instrument. The customer should acquire appropriate training documentation and fax the proof of training, when completed, to BioGenex Customer Support.

Recommendations / Customer comments:

Final Notes:

- The person assigned to training must be free from interruptions during the on-site training. It is preferable to train a “key” operator, and then let them train others as necessary.

- Inform the customer that they must be available on the days of installation. Complete installation of the instrument, calibration and training may take at least two days.

Recommendations / Customer comments:

Signatures

CUSTOMER:

BIOGENEX REPRESENTATIVE:

Date: _____

Date: _____

Reviewed by BioGenex Customer Support: _____

Reviewer’s Signature

Date

Appendix-8

Instrument Installation Instructions for i6000™ Diagnostics/Infinity Automated Staining Systems

1. Purpose

1.1 To define the process for installation of i6000™ Diagnostics/Infinity systems quickly and efficiently.

2. Scope

2.1 This procedure applies to installation of the i6000™ Diagnostics/Infinity systems.

3. Definition: N/A

4. Responsibility

4.1 Field Service Engineers, Application Specialists or appropriately trained BioGenex personnel have to install the unit.

4.2 Installers have to fill out the Installation Checklist Document No. 951-6124.1 in duplicate.

4.2.1 Provide copy of the completed checklist to the customer.

4.2.2 Provide original of the completed checklist to BioGenex Customer Support.

5. Related Documents:

5.1 Packing List for Accessories Box - i6000™ Diagnostics/Infinity systems. (Doc# 951-3323.1)

5.2 Instrument Installation Checklist for i6000™ (Doc# 951-6124.1)

5.3 Operator's Manual i6000™ Diagnostics System Doc No.953-0048.0

5.4 Operator's Manual i6000™ Research System Doc No.953-0049.0

6. Procedure:

6.1 Verify the received accessories against the Packing Lists. Contact Customer Support in case of short shipment and await instructions.

6.2 Remove the i6000™ Diagnostics/Research system from its crate and place it in the Lab on a flat surface that has sturdy support.

6.3 Ensure that the instrument is located with 6" clearance from the rear wall for ventilation.

6.4 Open the housing cover and remove the packing material inside the instrument.

6.5 Remove the M3x16 SHCS screw fixed to the left of the Z Head cover. This screw is fixed at BioGenex to ensure that the pipette adapter is secured firmly during the shipment.

6.6 Unpack the Computer and the monitor.

6.7 Ensure that the 110/230 Voltage Selector Switch located on the rear side of the computer is set per the local Mains Power System.

6.8 Connect the video cable to the Computer and the Monitor.

6.9 Connect the Keyboard and Mouse to the Computer.

6.10 Connect the Ethernet Cable (4500-00132; Cable Ethernet Cross Over) from the Computer to the port marked "Ethernet" on the connector panel of i6000™.

6.11 Connect the 9-pin barcode cable to the COM port on the Computer and the port marked "Scanner" on the connector panel of i6000™ (Applicable to i6000™ Diagnostics Systems only).

6.12 Install the Slide Racks, Reagent Racks and Pipette tips.

6.13 Install the Pipette Tip Waste Tray.

- 6.14 Fill the DI Water and Buffer carboys with at least 2L of DI Water.
- 6.15 Attach the carboy caps to their respective carboys and connect them to the liquid panel per the color coding (Buffer: White; DI Water: Green)
- 6.16 Attach drain connector (Elbow Fitting) to one end of the large drain tubing. Plug the elbow fitting into the connector on the liquid panel and the other end of the drain tubing to the Waste Carboy Cap and secure with a cable tie.
- 6.17 Place 4 slides in positions 1, 12, 49 and 60 on the slide racks and level the instrument. Adjust the leveling screws as required with the help of Level (6520-04116) available in the accessories kit.
- 6.18 Connect the power cords of i6000™, Computer, Monitor, Printer* and i500* Labeling System to a UPS System* of appropriate rating.

* Wherever supplied.

NOTE: Perform the setup of Inkjet Printer and i500 Plus - LIS Enabled Barcode Printer wherever supplied with the i6000™ system.

- 6.19 Turn ON the Computer and Monitor.
 - 6.20 Select the user "Autostainer". Type in password "optimax".
 - 6.21 Check the date and time and make changes if necessary.
 - 6.22 Turn ON the i6000™ Diagnostics/Research system. Ensure that the Emergency Stop switch is not pressed.
 - 6.23 Perform the following calibrations per procedure explained in Section 11 of the Operator Manual 953-0048.0 and 953-0049.0
 - 6.23.1 Multipoint Calibration
 - 6.23.2 Barcode Reader (Slide and Reagent).
 - 6.23.3 Vial Bottom
 - 6.24 Prime the Buffer and DI Water lines to verify the liquid flow.
 - 6.25 Check the drain flow and adjust the position of the Waste Carboy or the length of the drain tubing to enhance the flow.
 - 6.26 Setup the Inkjet Printer per the following procedure.
 - 6.26.1 Switch ON the printer.
 - 6.26.2 Install the printer cartridges.
 - 6.26.3 Connect the printer cable to the USB port of the Computer.
 - 6.26.4 Insert the CD provided by the Printer Manufacturer in the CD-ROM of the computer to install the printer drivers.
 - 6.26.5 Print a test page to ensure that the printer is ready.
 - 6.27 Perform a 4 Corner plus 2 Test per procedure explained in Appendix 5 of the Operator Manual 953-0048.0
 - 6.28 Drain out the left over DI water in the carboys, if any.
 - 6.29 Install the i500 Plus - LIS Enabled Barcode Label Printer per the following procedure.
 - 6.29.1 Connect the i500 Plus - LIS Enabled Barcode Label Printer to the USB port of the computer.
 - 6.29.2 Connect the Power adapter of the i500 Plus - LIS Enabled Barcode Label Printer to the UPS System*
- * Wherever supplied.
- 6.29.3 Insert the i500 Plus - LIS Enabled Barcode Label Printer CD in the CD-ROM of the computer to install the setup files.
 - 6.29.4 Install the printer drivers available in the i500 Plus - LIS Enabled Barcode Label Printer CD
 - 6.29.5 Restart the computer after the installation of drivers is complete.
 - 6.29.6 Click "Barcode Labeler" icon on the desktop.

- 6.29.7 Select i6000™ Automated Stainer from the list of products.
- 6.29.8 Click “Sync” to synchronize the Labeler with i6000™ Diagnostics/Research system. Ensure that the i6000™ system is turned ON while performing synchronization.
- 6.30 Train the customer’s key operator on the staining and diagnostic functions of the instrument.
- 6.31 Peel off the white backing from the Lamination Holder and stick it on Fab. Side Housing next to the serial number label.
- 6.32 Insert the laminated Daily/Maintenance Scheduler in the Lamination Holder from the top.
- 6.33 Peel off the backing from the Business Cards Holder and stick it on Fab. Side Housing next to the Lamination Holder.
- 6.34 Insert the Business Cards of the FSE/AS/Sales in the Business Cards Holder.
- 6.35 Fill out the Installation Checklist document 951-6124.1 and file copies of Checklist as mentioned in Section 4.2.